
FACELIFT & TONE

the beauty specialists

ENGLISH	4
FRANÇAIS	10
DEUTSCH	16
NEDERLANDS	22
ESPAÑOL	28
PORTUGUÊS	34
ITALIANO	40

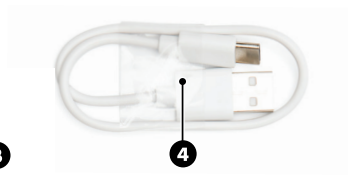WARNINGS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Do not place or store in direct sunlight. Store in a cool dry place.
- Keep out of the reach of children.
- If you suffer from a skin condition, consult your dermatologist before use
- For hygiene reasons we recommend for single person use only.
- Do not use this appliance on loose skinned parts of your body or any areas showing signs of: varicose veins, skin diseases (for example eczema, psoriasis, dermatitis), cuts and abrasions, irritation, warts, moles, skin tags, skin rash, hypersensitivity or sunburn, phlebitis, recent scar tissue, other skin conditions. We recommend you consult your physician before treatment if you suffer from any of the above conditions.
- Do not use on cuts or broken skin.
- Do not use on areas that have had cosmetic surgery.
- Do not leave the device unattended. Turn off before you leave the room.
- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Do not immerse product in water.
- Do NOT use this product if it is not working properly, or if it has been dropped or damaged.

KNOW YOUR UNIT

(see page 3)

- 1 Protective Cap
- 2 Treatment Heads
- 3 LED Lights (Blue, Green, Red)
- 4 USB Cable
- 5 Mini USB Socket
- 6 On/Off
- 7 Conductive Gel

INTRODUCTION

The FaceLift & Tone is a non-invasive at-home facial toner intended for cosmetic use and can be incorporated in to your daily/weekly beauty regime.

Your FaceLift & Tone works by delivering soft & gentle microcurrent pulses to the skin. The two contact points are designed to glide smoothly over the skin and to follow the natural contours of the face.

IMPORTANT - PATCH TEST

Before you start you must ALWAYS do a patch test with the conductive gel on an inconspicuous area of your skin for possible skin reaction. Wait 24 hours. If irritation occurs do not continue treatment.

CHARGING

- Plug the USB cable into the unit and then into a USB port.
- The LED Light will flash red to indicate charging, and turn green when fully charged. When fully charged the unit will automatically switch off after 30 minutes.
- It takes about 1.5 hours to fully charge.
- When fully charged the toner should operate for approximately two weeks of daily 20-minute treatments.

GETTING STARTED

1. Clean and dry the skin in the area to be treated.
2. Apply a generous amount of conductive gel to the Treatment Heads or to the skin in the treatment area.
3. Pressing the On/Off button for three seconds to turn on. The blue LED will indicate level 1 (Low) is selected. Press the On/Off button again to select Levels 2 and 3.

Treatment Levels:

- Level 1 – Low (Blue)
- Level 2 – Medium (Green)
- Level 3 – High (Red)

The Treatment Heads will warm up to 43°C.

4. Press the On/Off button for three seconds to turn off.

TREATMENTS

Having applied conductive gel slowly glide the Treatment Heads over the skin in the direction of the arrows (see treatment areas on page 4) with each movement lasting 10-15 seconds. Repeat 4-5 times.

- ➊ Cheeks - Run from the nasal lines in an upwards swooping direction over the cheek.
- ➋ Jaw Line - Follow the jaw line outward from the chin towards the ear.
- ➌ Under Chin - Glide upwards from the point where the chin meets the neck to the tip of the chin.
- ➍ Neck - Run from either side of the mid-line of the neck outwards. Repeat for the lower and upper sections of your neck.
- ➎ Forehead & Brow Line - Move upwards from the brow line to the hair line. Do this centrally and to the right and left sides of the forehead.

When used correctly the stimulation is soft and gentle - you will not necessarily feel the muscles contracting.

Important: Do not use for more than five minutes in any one area. **The unit will automatically turn off after 5 minutes.**

AFTER USE

- If you have very sensitive skin, you may notice redness after use. This should subside after 30 minutes.
- When not in use fit the Protective Cap.

CLEANING AND MAINTENANCE

- Before and after use, wipe the Treatment Head with a clean damp cloth, do not immerse the unit in water.

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly

SPECIFICATION

Model:
Input:
Battery

FALT

5V --- 500mA 2.5W

Li-Ion 3.7V 600mAh

Read these instructions carefully before use and then retain for future use.

This product is for indoor use only.

This product complies with current EU product legislation.

CUSTOMER CARE

If you have any questions, please contact us at:

online@riobeauty.com

GUARANTEE

To register your two year guarantee online, visit:

www.riobeauty.com/pages/warranty-registration

www.riobeauty.com/collections/consumables/products/conductive-gel-85ml

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© The Dezac Group Ltd 2023 Manufactured by The Dezac Group Ltd PO Box 17, Cheltenham, England GL53 7ET