
SCAN QR CODE

To access the full user manual

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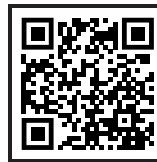
hairmax®
TRANSFORMATIVE HAIR REGROWTH

QUICK START GUIDE

LaserBand 82
comfortflex

QUICK START GUIDE

LaserBand 82 *comfortflex*

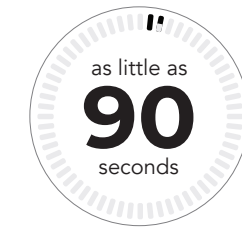as little as
90
seconds

3 times
per week


FOR MAXIMUM BENEFIT

- Scalp should be freshly washed (hair can be wet, dry or damp)
- Use your device prior to applying styling products (these can build up and impede the laser light)
- Regularly clean the laser area with a damp cloth. Do not use harsh detergents or cleaning products on any part of your device


CHARGING OPTIONS

Charge your LaserBand for 4-6 hours before initial use


Prior to charging your LaserBand, make sure your device is turned off.

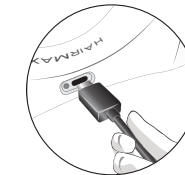OPTION 1 Plug the power cord into the charging stand and the other end into an electrical wall outlet. Place your device on the charging stand. The device indicator white light will flash during the charging process. A steady light indicates your device is fully charged.

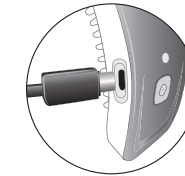OPTION 2 Place the power cord directly into the charging socket on the side of your device and the other end into an electrical wall outlet. The device indicator white light will flash during the charging process. Once charging is complete, the light will stop flashing and remain on.

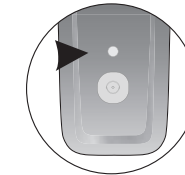The device indicator white light will flash during the charging process. A steady light indicates your device is fully charged.

Do not use your LaserBand while it is plugged into the power cord.

Prior to treatment, test your skin for sensitivity by shining the light on the inside of your forearm for 2 minutes. If you see any reaction, do not use your LaserBand and contact customer service.

TREATMENT METHOD

Your LaserBand will be used in 3 positions on your scalp. Leave it in place for a minimum of 30 seconds for each position. It will vibrate every 30 seconds to let you know when to move it to the next position.

To turn on your device, press the power button once to activate the red lasers.

- To start your treatment, press the power button once. Place your device at your hairline and leave it in place until it vibrates. This is position 1.
When you feel the vibration, slide your device backwards about 3 inches to position 2, the middle portion of your scalp.
- Use light pressure so the teeth can separate your hair. Leave it in position 2 until you feel the vibration again.
- Repeat the 30 second procedure on position 3. You have now completed today's treatment. Store your device in the cradle until your next treatment.

After treatment, remove your LaserBand from your head and turn off your device by pressing the power button once. Place your device in the charging cradle/stand and store for future use.

NOTE: The LaserBand 82 will turn off automatically after 10 minutes.

The LaserBand 82 is manufactured using precision components and should last for many years of normal use.

For detailed information about troubleshooting, proper use, warnings and cautions, access the full user manual by scanning the QR code located on the back panel. Please read the entire user manual before operating.

Your Hairmax device has a limited warranty against defects in materials or workmanship from Lexington International, LLC. Items that are damaged or abused are not covered under the warranty and will be charged at their regular repair rates. Any attempt to self-service your device, unless instructed by LEXINGTON INTERNATIONAL CUSTOMER SERVICE, will void your warranty. Please refer to our website for complete warranty terms.


To validate your warranty, register online at
www.hairmax.com/warranty

Laser Light - Avoid Direct Eye Exposure

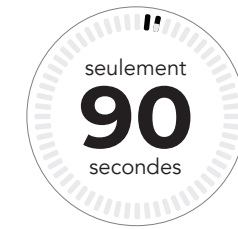Hairmax® laser devices are indicated to treat Androgenetic Alopecia, and promote hair growth in males who have Norwood-Hamilton Classifications of IIIa to V and in females who have Ludwig (Savin) I-4, II-1, II-2, or frontal patterns of hair loss and who both have Fitzpatrick Skin Types I to IV.

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