GUIDA RAPIDADual Action **18** LaserComb

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8
minutes
3 times
per week

FOR MAXIMUM BENEFIT


- Scalp should be freshly washed (hair can be wet, dry or damp)
- Use your device prior to applying styling products (these can build up and impede the laser light)
- Regularly clean your device with a slightly moistened cloth. Do not use harsh detergents or cleaning products on any part of your device

CHARGING OPTIONS

Charge your LaserComb for 2-3 hours before initial use

OPTION 1 Charge your device using the charging cradle and power cord. Prior to placing your device in the cradle, make sure it is turned off and the teeth face away from you. Plug the power cord into the cradle and the other end into an electrical wall outlet.

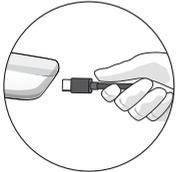OPTION 2 Charge your device by placing the power cord directly into the device and the other end into an electrical wall outlet.

The indicator light on the power button will flash during the charging process. Once charging is complete, the light will stop flashing and remain on.

When the battery is fully charged, the indicator light will be on. As the battery power diminishes, the light will blink slowly. When the battery is getting low and needs to be recharged, the light will start to flash rapidly. When this happens, charge your device.

If your device needs to be recharged during your treatment, allow it to fully charge before continuing.

Do not use your LaserComb while it is plugged into the power cord.

During use, the indicator light will indicate the battery level. A steady, non-flashing light indicates the battery is full.

Prior to treatment, test your skin for sensitivity by shining the light on the inside of your forearm for 2 minutes. If you see any reaction, do not use your device and contact customer service.

TREATMENT METHOD

Turn on your device by pressing the power button (lower button) until the power button light turns on. Next, press the laser button (upper button) until it lights up. That means your device is ready for use.

To start your treatment, place your device flat on your scalp so that both rows of teeth are touching the scalp.

Start at your hairline. Move your device from front to back of scalp in a slow combing motion, pausing every 1/2 inch (1 cm) for 4 seconds. When your device is on, it will beep and vibrate every four (4) seconds. This means that it's time to move the your device to the next area of the scalp. Do not lift your device off your scalp while combing from front to back. If you wish to turn off the beep/vibration during the treatment, press the top laser button again.

Continue in this manner from front to back and then side to side until you cover your entire scalp with multiple passes. Repeat these passes for the total treatment time of 8 minutes.

After treatment, turn off your device by pressing the power button (lower button) until the indicator light turns off. Store your device in the charging cradle for future use. You may leave the power cord plugged into the charging cradle at all times.

NOTE: The Actonis Dual Action 18 LaserComb will turn off automatically after 8 minutes.

The Actonis Dual Action 18 LaserComb is manufactured using precision components and should last for many years of normal use.

For detailed information about troubleshooting, proper use, warnings and cautions, access the full user manual by scanning the QR code located on the back panel. Please read the entire user manual before operating.

Your Actonis device has a limited warranty against defects in materials or workmanship from Lexington International, LLC. Items that are damaged or abused are not covered under the warranty and will be charged at their regular repair rates. Any attempt to self-service your device, unless instructed by LEXINGTON INTERNATIONAL CUSTOMER SERVICE, will void your warranty. Please refer to our website for complete warranty terms.

To validate your warranty, register online at
www.actonis.com/warranty

Laser Light - Avoid Direct Eye Exposure

Actonis™ laser devices are indicated to promote hair growth in females with Androgenetic alopecia who have Ludwig (Savin) Classifications I - II or frontal patterns of hair loss and who have Fitzpatrick Skin Types I - IV.

Actonis[™]
by hairmax®

SCAN QR CODE

To access the full
user manual

Or visit: www.actonis.com/usermanual

You can also call customer service free of charge to request a hard copy

CUSTOMER SERVICE

Direct line +1 561-417-0200

United States and Canada (toll-free) 1-866-527-3726

Authorized Representative in Europe (Regulatory Affairs Only)

EC	REP	EMERGO EUROPE Westervoortsedijk 60 6827 AT Arnhem The Netherlands
----	-----	--

Tel: +31 (0) 70 345-8570 Fax: +31 70 346-7299

Australia Sponsor

AU	REP	EMERGO AUSTRALIA 201 Sussex Street Darling Park, Tower II, Level 20 Sydney NSW 2000 Tel: +61 2 9006 1662
----	-----	--

Tel: +61 2 9006 1662

Brazil Sponsor

BR	REP	EMERGO BRAZIL Avenida Francisco Matarazzo 1.752, Salas 502/503 Água Branca São Paulo - SP CEP - 05001-200 Tel: +55 11 4800 9226
----	-----	---

Tel: +55 11 4800 9226

United Kingdom Sponsor

UK	RP	EMERGO UNITED KINGDOM Emergo Consulting (UK) Limited c/o Cr360 - UL International Compass House, Vision Park Histon Cambridge CB24 9BZ Tel: +44 (0) 1223 214 044
----	----	--

Tel: +44 (0) 1223 214 044

LEXINGTON INTERNATIONAL, LLC
1040 Holland Drive | Boca Raton, Florida | 33487 | USA
www.actonis.com

©2024 Lexington Intl., LLC. All Rights Reserved

LQH018DAIDXXX 00 0524