# AILORIA

## SHINE BRIGHT

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**EN** INSTRUCTION MANUAL

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
Version 2.0 - 10/2021

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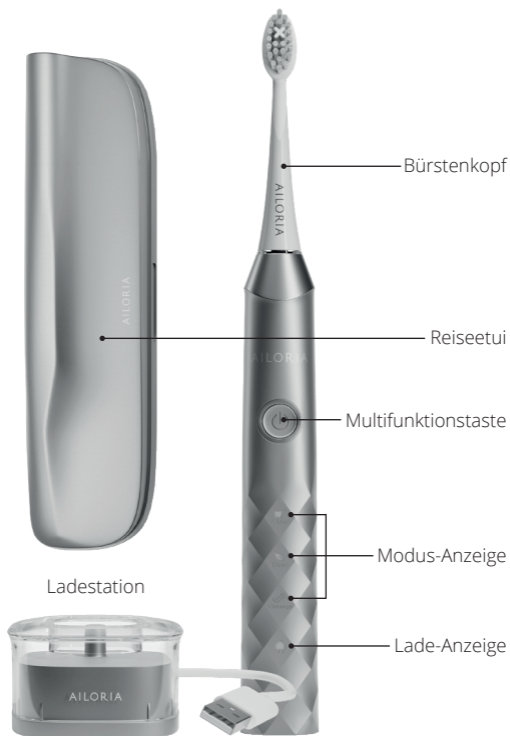
© VEHNS GROUP GmbH
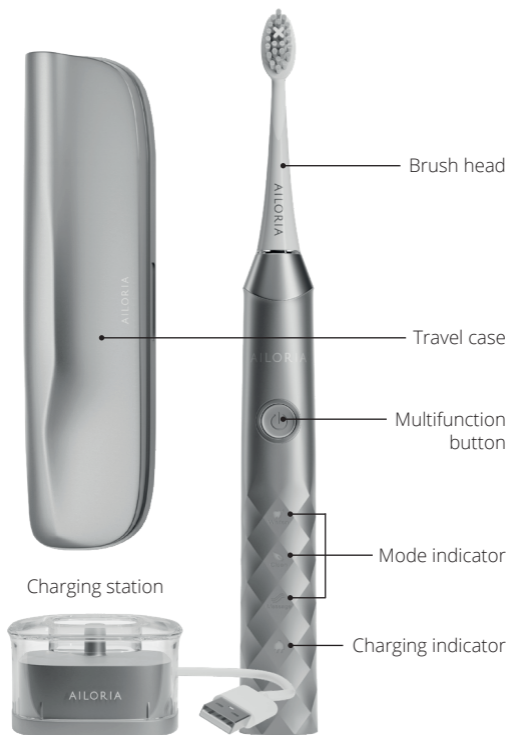
## **PRODUCT FEATURES**

The efficient sonic technology of the SHINE BRIGHT sonic toothbrush removes plaque and discolouration. With 31,000 brush head movements per minute and a wide bristle oscillation, the SHINE BRIGHT sonic toothbrush cleans effectively and thoroughly.

- 3 cleaning modes: Whitening, Clean, Massage
- 2 minute timer and 30 second quadpacer
- Replacement brush head with wavy bristles for gentle and effective care of teeth and gums
- Lithium battery and inductive charger with low power monitoring and overcharge protection
- Charging of the lithium-ion battery via USB
- Up to 21 days battery life with 2x daily use for 2 minutes each time
- Waterproof according to IPX7 standard

## **SAFETY INSTRUCTIONS**

- Use SHINE BRIGHT only for its intended use as described in this manual. This product is designed to clean teeth and gums. If you experience any pain or discomfort, please stop using SHINE BRIGHT and consult your dentist or doctor as this may be a sign of other health

problems.

- Consult your doctor before using SHINE BRIGHT if you have any health problems.
- If you have had any surgery on your teeth or gums during the last 2 months, please consult your dentist before using the sonic toothbrush.
- If excessive gum bleeding or pain occurs and persists during or after the use of SHINE BRIGHT, please consult your dentist.
- SHINE BRIGHT should not affect pacemakers in any way. However, for maximum precautions, we recommend that you consult your doctor before use.
- SHINE BRIGHT is only suitable for cleaning the teeth, gums and tongue and should not be used for any other purpose.
- Children under 8 years of age and persons with reduced physical or mental capabilities or lack of experience and knowledge can only use SHINE BRIGHT under adult supervision and instruction. However, cleaning and care of the product must only be carried out by adults with sufficient skills. Always keep the sonic toothbrush out of the reach of children.
- Do not use damaged brush heads at any time. Deformed bristles may break off during brushing. Replace the brush head every 3 months or sooner if signs of wear appear. Be aware that improper cleaning can damage the

bristles.

- Do not use the sonic toothbrush if it is damaged (brush head, handle, charger or USB cable). The sonic toothbrush and its individual components can only be repaired by the manufacturer.
- Disconnect the plug from the socket before cleaning the charging station.
- Never clean the sonic toothbrush, the brush attachments, the charger and the USB charging plug in the dishwasher.
- Only clean SHINE BRIGHT under running water. The charger must only be wiped with a damp cloth and not held under water.
- Do not use any attachments other than those recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or other bicarbonate, clean your brush head and rinse it thoroughly with water after each use. This will prevent plastic breakage due to prolonged exposure to these chemicals.
- Always keep the charger including the USB cable away from water as well as other liquids.
- Check that the specified voltage matches the local voltage before connecting the unit.
- The USB cable is permanently attached to the charging station and cannot be replaced. If it is damaged, please contact your contractor. Only have the charging station replaced with an

original spare part.

- Do not use the charger on hot surfaces.

## **COMMISSIONING**

Charge the sonic toothbrush before first use. Charging time is approximately 8 hours to charge the toothbrush before first use or fully discharged battery.

### **Charging**


- Connect the charging base (included) to a power source using a USB power adapter. The USB power supply unit is included in the scope of delivery depending on the selected version. If it is not included, any commercially available USB power supply unit can be used (5V 1A). Alternatively, SHINE BRIGHT can also be charged at other USB interfaces (e.g. notebook, PC, etc.).
- The charging process starts automatically as soon as the sonic toothbrush is placed on the charging station if it is connected to a power source.
- If the battery is only slightly charged, the charging plug symbol on the sonic toothbrush flashes ten times after switching off to remind you to charge the toothbrush.

- During the charging process, the charging plug symbol lights up red. Once charging is complete and the battery is fully charged, the light goes out.
- The toothbrush cannot be used during the charging process. If the toothbrush is placed on the charging station during operation, it switches off and goes into charging mode.
- When charging is complete, unplug the charger from the wall socket and store it in a dry place. Do not immerse the charger in water for cleaning, but wipe it with a damp cloth.

**Note:** It takes approx. 8 hours to fully charge the battery. The useful life of a full battery charge is 21 days with two minutes of use twice a day. If the red indicator light flashes or the product does not switch on, the battery power is insufficient and the product needs to be recharged.

## **Installing the brush head**

Push the brush head onto the metal shaft of the toothbrush until it is secure. To remove, pull the brush head firmly upwards.

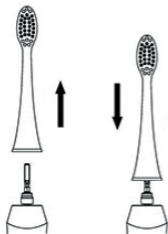
Dentists recommend replacing the brush heads every 3 months.

## APPLICATION

- Moisten the bristles of the sonic toothbrush, apply toothpaste and switch on the toothbrush by pressing the on/off switch when the brush head is in your mouth.
- Throughout brushing, move the brush head at a slight angle (approx.  $45^\circ$ ) and slowly over the teeth and chewing ridges without applying pressure. The brush head should be placed with one half on the gums and the other half on the teeth to ensure cleaning also on the important transition between gums and teeth. The high-quality wavy bristles, which are adapted to the shape of the teeth, together with the sonic movements clean your teeth in a gentle way.

- Brush your teeth for 30 seconds on each of the quadrants. The interval timer indicates when the quadrant should be changed by briefly interrupting. The smart timer ensures that the two-minute brushing time recommended by the dentist is adhered to.

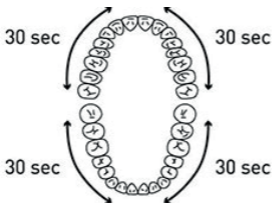- When you switch SHINE BRIGHT on for the first time with the On/Off button, your toothbrush automatically starts in Whitening mode.
- Press the power button and switch on the toothbrush. Press the On/Off button again within 3 seconds to switch to the next mode. Repeated pressing activates all 3 modes continuously. After the third and last mode, the toothbrush switches off.
- The indicator lights on the toothbrush show you which mode is activated.
- To switch off the sonic toothbrush, press the On/Off button 3 seconds after switching it on. If you press the button too soon, you will switch to another brushing mode.

- Due to the integrated memory function, the toothbrush will automatically start in the same mode you left off with the next use.

## CLEANING MODES

The SHINE BRIGHT sonic toothbrush has 3 cleaning modes:

**Whiten:** This mode provides intensive cleaning of the teeth, removing plaque and whitening them with regular use. This mode is suitable for users with insensitive teeth.

**Clean:** This mode is for daily thorough cleaning. The vibrations of the bristles are slightly weaker than in Whiten mode. This mode is also suitable for users with sensitive teeth.

**Massage:** The vibration of the bristles and the intensity vary to promote blood circulation to the gums. Use this mode to polish your teeth and massage your gums.

## CLEANING & STORAGE

- Remove the brush head and rinse the metal shaft with warm water. Do not press against the rubber seal on the metal shaft with sharp

objects as this may damage it.

- Wipe the entire surface of the handle with a damp cloth.
- Rinse the brush head and bristles under running water after each use. Thanks to the IPX7 protection standard, SHINE BRIGHT (without the charger) is waterproof and can also be used in the shower.
- Remove the brush head from the handpiece once a week and rinse it thoroughly with warm water.
- To clean the charger, remove the transparent cover of the charger and clean it with a damp cloth. Never immerse the charger in water or rinse it under running water.
- For optimal cleaning results, you should replace the brush heads every 3 months when signs of wear are evident. Please only use replacement brushes recommended by the manufacturer.
- If you do not intend to use the SHINE BRIGHT for a longer period of time, please disconnect the USB cable or the plug from the socket and clean the toothbrush. Please store it in a cool, dry place out of direct sunlight.

## **TROUBLESHOOTING**

### ***The toothbrush will not switch on***

SHINE BRIGHT may need to be charged.

### ***Tingling sensation during use***

Due to the sound vibrations, you may experience a tingling sensation when using SHINE BRIGHT for the first time(s). This feeling will decrease over time as you get used to the vibrations.

### ***The sonic vibrations seem weak***

The bristles may be pressed too hard against your teeth. Gently place the brush head on your teeth.

## **ACCESSORIES & REPLACEMENT BRUSHES**

Replacement brush heads, charging station and charging plug are available from your contract partner or directly from the manufacturer: [www.ailoria.de](http://www.ailoria.de)

## **WARRANTY & SERVICE**

Defects due to material and/or workmanship will be repaired or replaced by AILORIA under warranty if proof of purchase is obtained during the warranty period. Contact the Customer Service Centre at [service@vehnsgroup.com](mailto:service@vehnsgroup.com).

Items not covered by the warranty:

- Brush heads
- Damage caused by misuse, abuse, neglect, alteration or unauthorised repair.
- Normal wear and tear including scratches, discolouration or fading.

## **TECHNICAL DATA**

Toothbrush power: 3 W

Charger power: 1.5 W

Toothbrush voltage: DC 3.7 V

Charger voltage: DC 5 V

Battery capacity: 700 mAh

## **DISPOSAL**

Do not dispose of the product with normal household waste. Follow local regulations for separate disposal of electrical and electronic products.

Proper disposal will help prevent any negative impact on the environment.

## **MANUFACTURER**

VEHNS GROUP GmbH  
Theatinerstrasse 40-42  
80333 Munich, Germany

service@vehnsgroup.com  
www.vehnsgroup.com www.ailoria.de


MADE IN CHINA

Subject to alterations and errors

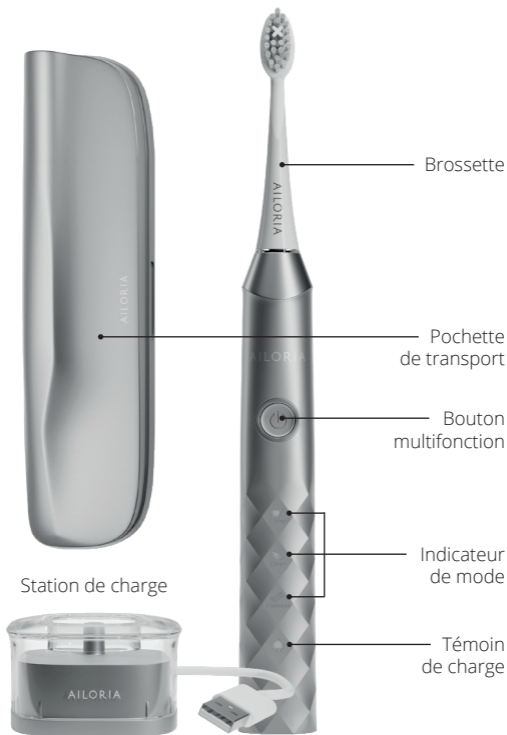The manufacturer accepts no liability for damage caused by incorrect, improper or incompatible use of the product.

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VEHNS GROUP GmbH, Theatinerstrasse 40-42