ORDO

SONIC+
ELECTRIC TOOTHBRUSH

Important Information

Wichtige Informationen

Une information important

Información importante

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Önemli bilgi

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Viktig information

Vigtig information

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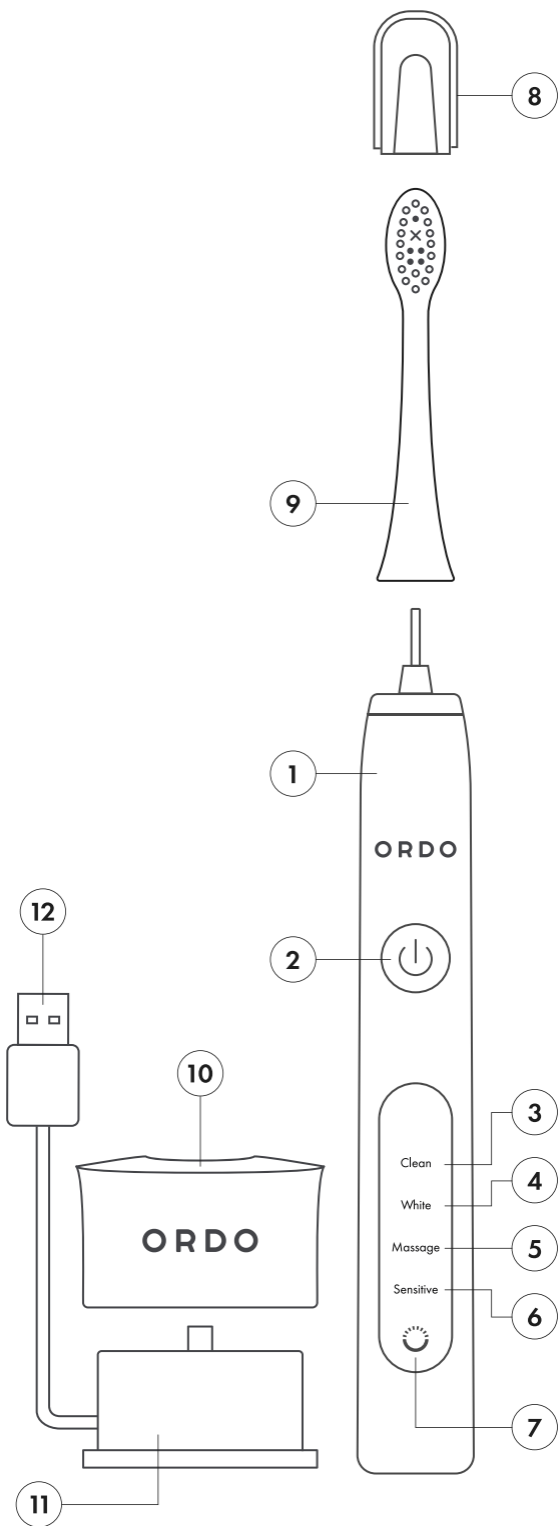
ordolife.com

Ordo, Manchester, WA14 4DR, UK.

Authorised Rep Compliance Ltd., Ground Floor,
71 Lower Baggot Street, Dublin, D02 P593, Ireland.

English

Technical Data

Model number: SP2000-CG, SP2000-RG, SP2000-WS, SP2000-PV, SP2000-MG, SP2000-IB

Charging Base

Input: DC 5V, 1.5W

Protection class: III

IP code: IPX7

Dimensions (LxWxH): 6.2cm x 6.2cm x 4.3cm

Toothbrush

Rechargeable battery: 3.7V, 700mAh

Operating time: 56 cycles/full charge

Timer: 0.5 second pause after 30 seconds, deactivation after 2 minutes.

Dimensions (LxWxH): 2.6 x 2.7 x 25.4cm

Warranty

This product is covered by a 2 year warranty. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty, please be sure to contact our customer service team for instructions. For more information please visit ordolife.com/warranty

60-Day Money-Back Guarantee

Satisfaction guaranteed – if you don't love your brush within 60 days we'll give you a refund. For more information on our guarantee please visit ordolife.com/guarantee

Customer Service

For more information on Ordo's products and services please visit ordolife.com

If the device is damaged, defective, in need of repair or you require assistance please contact our customer service team:

Email: hello@ordolife.com

Support: help.ordolife.com

Our team are always happy to help.

Recycling & Disposal

At Ordo we are always looking for solutions to reduce packaging and product materials that cannot be recycled effectively. We ask every member of the Ordo family to help us by recycling responsibly.

For more information, please visit ordolife.com/recycling

Do not dispose of the product with normal household waste at the end of its life. Take it to an official collection point for recycling and proper disposal.

This product contains a built-in, lithium-ion rechargeable battery. The battery must only be removed by a qualified professional when the appliance is disposed of.

Your Sonic+ Toothbrush

The list below refers to the image on page 2.

1. Handle
2. Power/Mode button
3. Clean mode indicator
4. White mode indicator
5. Massage mode indicator
6. Sensitive mode indicator
7. Battery indicator
8. Brush head travel cap

9. Brush head
10. Charging base top cover
11. Charging base dock
12. USB-A connector

Using Your Sonic+ Toothbrush

If you're using your Sonic+ toothbrush for the first time, it's normal to feel more vibration than if you were using a manual or battery-powered toothbrush. Some people notice a slight tickling or buzzing sensation – this will pass in time. We recommend using the Sensitive mode for the first 1-2 weeks, then trying the Clean mode once you've become accustomed to the sensation.

Your Sonic+ toothbrush has four modes: Clean, White, Massage and Sensitive. Each mode can be selected by pressing the power button repeatedly until you reach your desired mode. Here is a breakdown of each mode:

Clean: perfect for powerful everyday cleaning.

White: high power mode for removing stubborn surface stains.

Massage: pulses of vibration stimulate blood flow to the gums.

Sensitive: ideal for those with sensitive teeth or gums.

Your Sonic+ remembers the last mode you used, so next time you use it, all you need to do is press the power button and you're good to go.

It's common for first-time users to apply too much pressure. The benefit of an electric toothbrush is that it does the brushing for you, so you only need to apply gentle pressure. Please follow the instructions below for the best brushing experience.

This toothbrush is safe to use on braces, dental restorations (fillings, crowns, bridges, implants), and cosmetic dental work (veneers, composite bonding), as long as the dental work is properly adhered and not compromised. Please ask your dentist for advice if you're unsure in any way.

Brushing instructions:

1. Wet the bristles and apply a pea-sized amount of toothpaste.
2. Place the toothbrush bristles at a 45-degree angle to the teeth, pressing just enough to make the bristles reach the gum line.
3. Press the power/mode button to switch on the toothbrush and select your desired mode.
4. Remember not to press too hard when brushing, let the vibrating bristles do the work of clearing away plaque and food debris – no need to scrub as you might with a manual toothbrush.
5. Guide the bristles across and between your teeth, using small back and forth motions.
6. To ensure you thoroughly brush all surfaces of your teeth, split your mouth into four quadrants: top teeth outer surfaces, top teeth inner surfaces, bottom teeth outer surfaces, bottom teeth inner surfaces. Use the inbuilt Quad-Pacer to guide you between each of the quadrants – every 30 seconds the vibrations will pause briefly to let you know it's time to move onto the next quadrant.
7. After 2 minutes the toothbrush will automatically switch off, letting you know that you've brushed your teeth for the dentist recommended time. Of course, you can switch the toothbrush back on to brush for longer if you'd like.
8. If at any time before the end of a 2-minute cycle you want to turn the toothbrush off, simply press the power/mode button.
9. Once you're done with brushing, spit out any excess toothpaste, but do not rinse your mouth with water. This allows the ingredients of the toothpaste to better protect your teeth.
10. Thoroughly rinse the bristles, brush head and handle, ensuring you wash away any toothpaste residue.

Charging Your Sonic+ Toothbrush

1. When the rechargeable battery is low, the battery indicator

will flash red when you finish a brushing cycle or power off the toothbrush.

2. Ensure the entire toothbrush is fully dry, paying particular attention to the bristles and the area around the bottom of the toothbrush.
3. Connect the USB-A connector of the Charging Base to a power source, then place the toothbrush into place, ensuring the front of brush is facing the same direction as the front of the Charging Base (the side with the "ORDO" logo).
4. The battery indicator will light up red to show the toothbrush is being charged.
5. Charge the toothbrush for 8 hours or until the charging light turns off.
6. Switch off the power source and disconnect the Charging Base.

Troubleshooting

Problem	Possible Cause	Solution
The toothbrush isn't charging.	The charging base may be defective.	Try other sockets or USB ports and charge the toothbrush for at least 8 hours.
	The USB connector may be defective.	
	The charging base is defective.	Replace the charging base within the warranty period. Contact our customer service team.
The toothbrush operates with only a small amount of vibration.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours.
The toothbrush does not run for 2 minutes.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours and make sure that the toothbrush is properly placed on the charging base.

WARNINGS & CAUTIONS

Safety Instructions

The toothbrush is designed exclusively for cleaning teeth. It is only intended for private use and is not suitable for commercial purposes. Only use the toothbrush as described in these operating instructions. Any other use is deemed improper and may result in damage to property or even personal injury. The toothbrush is not a child's toy. Ordoclub Ltd accepts no liability for damage caused by improper or incorrect use.

Risk of electric shock

A faulty electrical installation or excessive mains voltage may result in an electric shock.

- Only connect the charging base if the mains voltage of the socket corresponds to the specification on the nameplate.
- The main cord of the charging base cannot be replaced. In the event of damage, the charging base must be disposed of.
- Only connect the charging base to an easily accessible socket so that you can quickly disconnect it from the mains in case of a failure.

- Do not operate the toothbrush and the charging base if they exhibit visible damage or the main cord or USB connector is defective.
- Do not open the housing; instead, have a qualified professional perform any repairs. Contact a qualified workshop for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection of the device or incorrect operation.
- Only parts that comply with the original device information may be used for repairs. Electrical and mechanical parts, which are essential for providing protection against sources of danger, are located in this toothbrush.
- The charging base must not be operated with an external timer or with a separate tele-control system.
- Do not immerse the charging base or the main cord in water or other liquids.
- Never touch the USB connector with damp hands.
- Never pull the USB connector out of the socket by the main cord.
- Never use the main cord as a carrying strap.
- Keep the charging base, toothbrush and USB connector away from open flames and hot surfaces.
- Do not kink the main cord and do not lay it over sharp edges.
- Only use the charging base in rooms located indoors.
- Never setup the charging base such that there is a risk of it falling in a bath tub or sink.
- Never attempt to grab the charging base if it has fallen into water. In such a case, immediately disconnect the USB connector.
- Only charge the toothbrush with the enclosed charging base.
- Make sure that children do not insert any objects into the charging base or the hand piece.
- Disconnect the USB connector when you are not using the charging base, after the charging process, before you clean the toothbrush or in the event of a fault.

Risk of chemical burns

Fluid leaked from the rechargeable battery may lead to chemical burns.

- Do not allow acid from the rechargeable battery to come into contact with skin, eyes or mucous membranes.
- In the event of contact with acid from the rechargeable battery, immediately flush the affected areas with plenty of clean water and seek medical help.

Danger of explosion

Improper handling of rechargeable batteries may cause an explosion.

- Never short circuit rechargeable batteries.
- Never dismantle rechargeable batteries.
- Never throw rechargeable batteries in fire.
- Never expose rechargeable batteries to excessive heat or direct sunlight.
- Do not in any case open the rechargeable batteries.

Danger for children, persons with impaired physical, sensory or mental capacities

e.g. partially disabled persons, older persons with reduced physical and mental capacities, and those lacking experience and knowledge (e.g. older children).

- This toothbrush may be used by children aged 8 years and above, persons with impaired physical, sensory or mental capacities, and those lacking experience and knowledge if they are supervised or have been instructed in how to safely use the toothbrush and have understood the risks associated.
- Cleaning and user maintenance must not be performed by children.
- Keep children under 8 years of age away from the toothbrush.

- Do not leave children unattended while using the toothbrush.
- Do not allow children to play with the toothbrush. Children could swallow small parts and choke.

Risk of injury

Handling the toothbrush improperly may result in injury.

- Lay the main cord so that it does not pose a tripping hazard.
- Immediately replace worn or damaged brush heads to prevent injury to the mouth and tooth area.

Risk of damage

Improper handling of the toothbrush may damage it.

- Do not drop the hand piece. If this does happen, do not use the toothbrush any more. Have it checked by an authorised service provider.
- Turn the toothbrush off when you change the brush head.
- Place the charging base on an easily accessible, level and dry surface.
- Ensure that the main cord, the charging base and the toothbrush do not come into contact with hot surfaces or open flames.
- Stop using the toothbrush if its plastic parts exhibit cracks or breaks or are deformed.

First Use

Checking the toothbrush and package contents.

Risk of damage

If you are not cautious when opening the packaging with a sharp knife or other pointed object, you may damage the toothbrush and injure yourself.

- Be very careful when opening the packaging.
1. Take the toothbrush and its accessories out of the packaging.
 2. Check to make sure that all the items are included.
 3. Check whether the toothbrush or the individual parts are damaged. If this is the case, do not use the toothbrush. Contact Ordo using the customer service contacts provided.

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Problem	Possible Cause	Solution
Lo spazzolino non si carica.	La base di ricarica è difettosa.	Sostituire la base di ricarica entro il periodo di garanzia. Contatta il nostro servizio clienti.
The toothbrush operates with only a small amount of vibration.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours.
The toothbrush does not run for 2 minutes.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours and make sure that the toothbrush is properly placed on the charging base.

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