


**ORDO**  
HYDRO SONIC  
WATER FLOSSER

<b>English</b>	<b>3</b>	<b>Svenska</b>	<b>31</b>	<b>Tiếng Việt</b>	<b>59</b>
<b>Deutsch</b>	<b>7</b>	<b>Dansk</b>	<b>35</b>		
<b>Français</b>	<b>11</b>	<b>Polski</b>	<b>39</b>		
<b>Español</b>	<b>15</b>	<b>Eesti keel</b>	<b>43</b>		
<b>Nederlands</b>	<b>19</b>	<b>Русский</b>	<b>47</b>		
<b>Italiano</b>	<b>23</b>	<b>51</b>	<b>العربية</b>		
<b>Norsk</b>	<b>27</b>	<b>汉语</b>	<b>55</b>		

1. Floss tip
2. Floss tip collar
3. On-demand button
4. Power button
5. Mode selection button
6. Main body
7. Water tank body
8. Water tank base

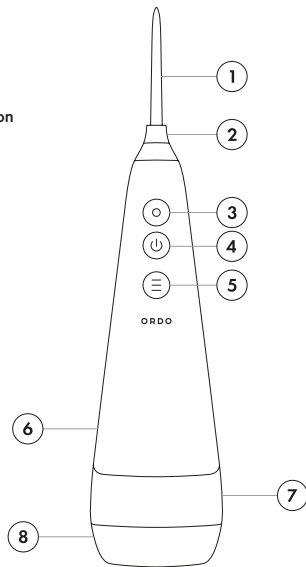
9. Mode indicator panel
10. Floss tip release button
11. USB-C charging port
12. Water tank release button
13. Water tank cap
14. Water tank base release button

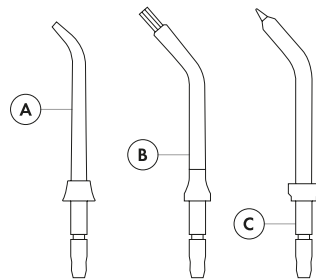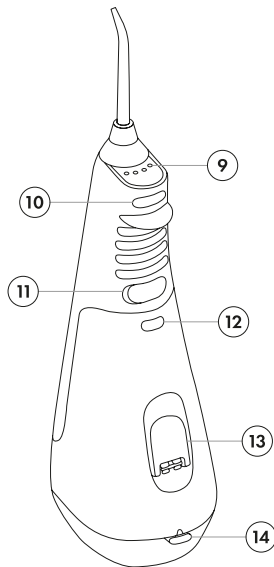
- A. **Premium Clean** – Ideal for excellent everyday cleaning. Suitable for most people.
- B. **Orthodontic** – Tapered bristles help clean around braces, bridges, implants and other dental work.
- C. **Periodontal** – Designed to rinse periodontal pockets, helping to protect against gum disease.


Orthodontic & Periodontal tips sold separately. Our full selection of floss tips are available at [ordolife.com](http://ordolife.com)

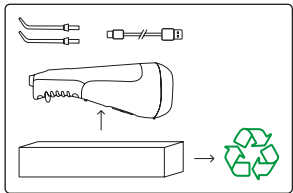**(1)** Remove all the contents from the box and check you have all parts specified. Make sure you recycle the inner and outer boxes after use.

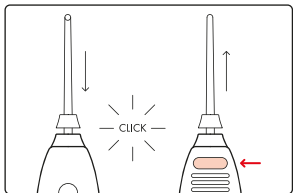
**(2)** Simply push the floss tip into the slot at the top of the water flosser, until it clicks into place. To eject, press the floss tip release button on the back of the water flosser, then pull the floss tip out.

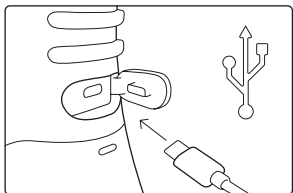
**(3)** Before using your water flosser, you should charge it for 6 hours. Open the USB-C port cover, then plug in the USB-C connector of the charging cable. Connect the USB-A connector to a power source. Ensure the water flosser is not on prior to charging.


**MAKE SURE THE USB-C PORT IS DRY BEFORE CHARGING!**

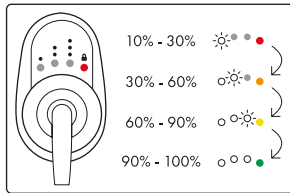
**(4)** During the charging cycle, the lights on the mode indicator panel will follow a sequence as shown in the diagram, to indicate the water flosser's charging status. When the coloured light turns green, the water flosser is fully charged. A full charge will last 30 days when using the water flosser for 2 minutes each day.

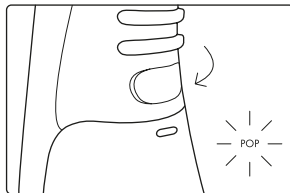
**(5)** Once you have charged the water flosser, remove the USB-C charging cable and securely close the cover. Make sure that this is closed properly to avoid water getting inside.

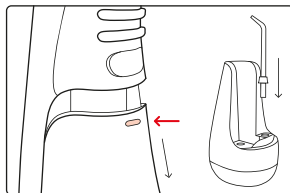**(6)** For travel convenience, you can store your Premium Clean floss tips inside the water tank, as shown in the diagram. Simply push the water tank release button at the top of the water tank and slide off. Inside you will find two slots to securely store your Premium Clean floss tips. To reattach the tank, slide back onto the main body until it clicks into place.

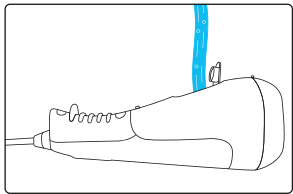
**(1)** To fill your water flosser, flip open the cap on the back of the water tank and fill with clean, lukewarm water (under 40°C). You could also add a few drops of Ordo Mouthwash Concentrate for extra freshness. If you use mouthwash, run a full tank of fresh water through the flosser after use, to avoid residue build-up.

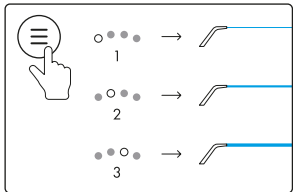
**(2)** Short press the mode selection button to scroll through the 3 pressure settings:

- 1 - Sensitive Mode (low) – recommended for new users.
- 2 - Daily Clean Mode (medium)
- 3 - Power Clean Mode (high)


You can change modes whilst using your water flosser.

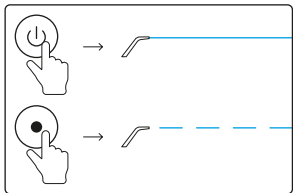
**(3)** Lean over the sink, place the floss tip in your mouth and part your lips slightly to allow water to flow out. Press the power button to begin flossing. Alternatively, press the on-demand button to control the flow of water – press and hold for a continuous flow and release to stop.

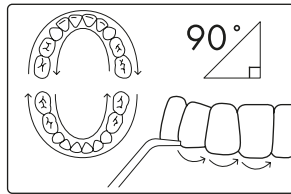
**(4)** Hold the floss tip at a 90° angle to your gum line and pause in between each tooth to ensure a thorough clean. Clean along the gum line as well as between teeth. Clean inside as well as outside surfaces. To get the most from the water tank, do not tilt the water flosser back during use. To help protect your water flosser, it will automatically turn off after running for 2 minutes.

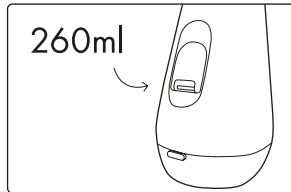
**(5)** Your water flosser has a 260ml tank capacity, which provides 55 seconds of use. We recommend 2 full tanks per flossing session, one for the top teeth and one for the bottom. Make sure the tank is fully emptied after use and open the cap to let it air dry. Wipe down the floss tip and main body.

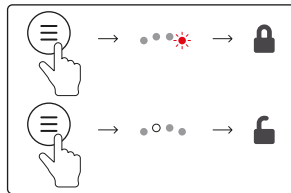
**(6)** To deactivate your water flosser for travel or storage, press and hold the mode selection button until the red light begins to flash, at which point the Lock Mode is activated. In Lock Mode the lock light will flash 5 times when you short press any button. Simply press and hold the mode selection button again to unlock the device – the last used mode will illuminate once your flosser is unlocked.

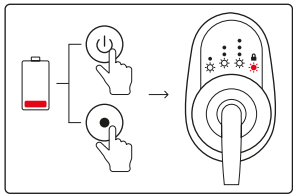
**(1)** At 30% battery, all the lights on the mode indicator panel will flash 5 times after use. When the water flosser has no power, the lights will also flash 5 times, but the device will not turn on – you need to recharge to a sufficient level before you can continue using the water flosser.

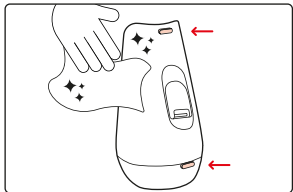
**(2)** To clean the water tank, first remove it from the main body by pressing the release button at the top of the tank and sliding off. Remove the base from the tank by pushing the release button at the bottom and pulling open. You can clean all parts of the tank with warm, soapy water or place in the dishwasher.

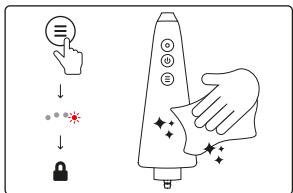**(3)** To clean the main body of the water flosser, ensure it is switched off and disconnected from any power source. We recommended activating the Lock Mode for cleaning. To avoid scratching, use a soft, dry or damp cloth to clean the main body.

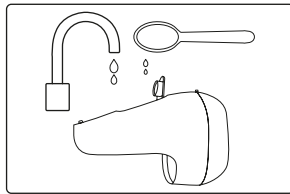**(4)** To clear hard water deposits from within your water flosser, fill the tank with water and add 1 tablespoon of white vinegar. Run the water flosser until the tank is empty, then rinse thoroughly with clean, warm water. Repeat this process every 1 to 3 months, depending on the hardness of your water.

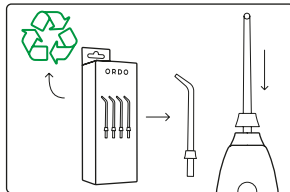**(5)** It is recommended that you replace floss tips when they become worn or dirty, this is usually every 6 months. Replacement tips as well as Orthodontic and Periodontal tips are available at [ordolife.com](http://ordolife.com)

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