(EN) USE & CARE GUIDE
FOR YOUR SAFETY AND CONTINUED ENJOYMENT OF THIS PRODUCT, ALWAYS READ INSTRUCTION LEAFLET CAREFULLY BEFORE USING.
The product is intended for private use only and is not intended for commercial use or for use in other applications. The product is designed for the neck and face area, do not use the product on other parts of the body.

SPECIFICATION:


Name: STYLPRO™ Neck & Face Smoother
Model: NS01A
Device Size: 89mm(L) x 47mm(W) x 120mm(H)
Materials: Coated ABS Plastic
Charging Voltage: 5.0V Nominal Voltage: 3.7V
Maximum Power: 2.5W

INCLUDED CONTENTS:

- 1 x STYLPRO™ Face and Neck Smoother
- 1 x Protective Cap
- 1 x Charging Cable
- 1 X Instruction Booklet

PARTS IDENTIFICATION:

- 1.** Message Head **2.** Display **3.** Power/Mode Button **4.** Red/Green/Blue Light Indicator **5.** battery charging indicator **6.** USB charging point.

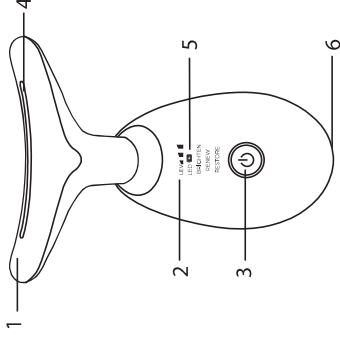CAUTION:

1. Use the device once daily.
2. We do not advise using this device if you: have undergone facial plastic surgery, have a metal implant, are pregnant, or if you are currently on your menstrual cycle.
3. Consult your medical professional prior to use if: You are epileptic, photophobic, or light sensitive, if you are pregnant or have had skin cancer in the last 12 months.
4. When the device needs cleaning, ensure that the power is off. Clean with a moist wipe, paper towel or cloth with a small amount of soap and wipe the head. Avoid contact with water and other metals. Do not use corrosive solutions such as alcohol, oils, or nail polish remover to wipe the device as this may lead to cracking, discolouration, or other damages to the device.
5. Store in a cool dry place. Do not store in a bathroom or room with humid air, which may damage the device.

CHARGING:

During the charging process, the display screen battery symbol flashes and the battery symbol stops flashing after charging is completed.


MAINTENANCE:

1. Please wipe off any excess skin care products which may accumulate in the massage head.
2. Wipe clean with a damp cloth only.
3. Do not submerge in water.
4. Please don't use cleaning product to clean.
5. Please store it in a cool dry place.
6. Please don't use if it is damaged.

DANGER: TO REDUCE THE RISK OF ELECTRIC SHOCK: ALWAYS TURN OFF BEFORE CLEANING OR STORING.

HOW IT WORKS:

Brighten mode
Brighten mode uses low levels of vibrations and **green** LED light therapy to help clear

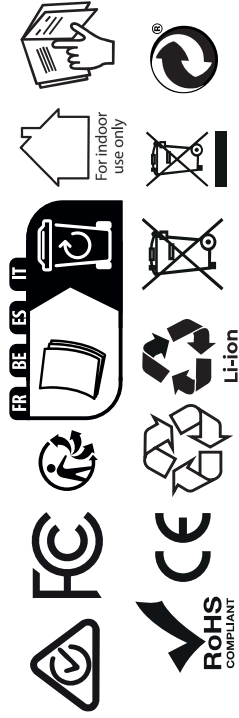©2025 STYLIDEAS Ltd. All rights reserved. STYLIDEAS owns trademarks, registered trade marks, patents, patents pending, design registrations and protection of the copyright of its products.
www.stylideas.com/lp. Distributed by STYLIDEAS Ltd,
UK: Amshold House, IG10 2RW, UK. +44(0)20 3225 5559
teamtom@stylideas.com.

EU: OBELIS S.A, Bd. Général Wahis 53, B-1030 Brussels, Belgium
AU: 6/7 Smeaton Grange Rd, Smeaton Grange, NSW.
USA: 2020 Ocean Ave Ronkonkoma NY, 11779

Made in China.

NS01A

the skin, working to reduce the appearance of hyperpigmentation and age spots on the face and neck.

Benefits:

- Promotes clearer skin.
- Helps reduce the appearance of hyperpigmentation and age spots.

Renew mode

Renew mode uses a medium level of vibrations and **red** LED light to encourage the production of collagen and elastin: naturally occurring proteins in the body that encourage firmer skin. This helps to reduce finer lines, wrinkles, and scarring.

Benefits:

- Encourages the production of collagen and elastin.
- Helps to make skin firmer, reducing the appearance of wrinkles, lines and scarring.

Restore mode

Restore mode uses a higher level of vibration, **blue** LED light and heat to help reduce acne, inflammation, and redness, preventing oil build up and reducing the occurrence of black heads in the skin.

Benefits:


- Helps reduce: Inflammation, Acne, Redness, Oil build-up, Black heads.

HOW TO USE:

Device movement

When using the device on your neck, we recommend these movements:

1. From the chin upwards towards your outer cheek
2. From the upper cheek upwards around your temples
3. From the bottom of your neck upwards towards your chin
4. From the bottom of your forehead towards your hairline

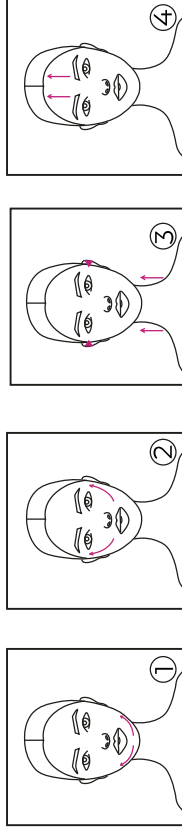Brighten mode

1. Hold down the power button for 1 second to turn on the device and enter the BRIGHTEN mode.
2. Apply moisturiser, serum, or lotion to the skin, massage until fully absorbed.
3. Ensure your fingers are touching the silver strip around the handle or the silver ring around the neck to enable the touch sensor.
4. The device will begin vibrating once the tip touches your skin.
5. Glide the massage head along your face and neck for at least 3 minutes. To target key areas with hyperpigmentation and age spots, we recommend holding the device on those areas for an extra 2-3 minutes.

Renew mode

1. Hold down the power button for 1 second to turn the device on, then quickly press the power button ONCE to enter the RENEW mode.
2. Apply moisturiser, serum, or lotion to the skin, massage until fully absorbed.
3. Ensure your fingers are touching the silver strip around the handle or the silver ring around the neck to enable the touch sensor.
4. The device will begin vibrating once the tip touches your skin.
5. Glide the massage head along your face and neck for at least 2 minutes. To target key areas with wrinkles and lines, we recommend holding the device on those areas for an extra 1 ½ minutes.

Restore mode

1. Hold down the power button for 1 second to turn the device on, then quickly press the power button TWICE to enter the RESTORE mode.
2. Apply moisturiser, serum, or lotion to the skin while you wait for the device to heat up. This will take at least 2 minutes. Massage until fully absorbed.
3. Ensure your fingers are touching the silver strip around the handle or the silver ring around the neck to enable the touch sensor.
4. The device will begin vibrating once the tip touches your skin.
5. Glide the massage head along your face and neck for at least 2 minutes. To target key areas with acne and inflammation, we recommend holding the device on those areas for an extra minute.

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(AU)
OUR GOODS COME WITH GUARANTEES THAT CANNOT BE EXCLUDED UNDER THE AUSTRALIAN CONSUMER LAW. YOU ARE ENTITLED TO A REPLACEMENT OR REFUND FOR A MAJOR FAILURE AND COMPENSATION FOR ANY OTHER REASONABLY FORSEEABLE LOSS OR DAMAGE. YOU ARE ALSO ENTITLED TO HAVE THE GOODS REPAIRED OR REPLACED IF THE GOODS FAIL TO BE OF ACCEPTABLE QUALITY AND THE FAILURE DOES NOT AMOUNT TO A MAJOR FAILURE.

(USA)
One (1) Year Limited Consumer Product Warranty
STYLIDEAS LTD ("STYLIDEAS") warrants that the product shall be free from defects in material and workmanship for a period of 1 year from the date of original purchase ("Limited Warranty Period"), provided the product was used in accordance with its use and care instructions (e.g., in the intended environment and under normal circumstances).

Your Limited Warranty Benefits
During the Limited Warranty Period and subject to this limited Warranty Policy, STYLIDEAS will, in its sole and exclusive discretion, either (i) refund the purchase price if the purchase was made directly from the online STYLIDEAS store, (ii) repair any defects in material or workmanship, (iii) replace the product with another product of equal or greater value, or (iv) provide store credit in the amount of the purchase price.

Claiming Your Limited Warranty Service in 5 Simple Steps:
1. Make sure your product is within the specified limited warranty period.
2. Make sure you have a copy of the invoice and order ID or proof-of-purchase.
3. Make sure you have your product. DO NOT dispose of your product before contacting us.
4. Contact our Customer Support Team via +44(0)20 3225 5559 teamtom@stylideas.com
5. Once our Customer Support Team has approved your request, please return the product with a copy of the invoice and order ID.

BATTERIES:
THE PRODUCT CONTAINS LITHIUM-ION BATTERY. RISK OF FIRE AND BURNS, DO NOT OPEN, CRUSH, HEAT ABOVE 60°C/140 °F OR INCINERATE. DO NOT SHORT CIRCUIT. IF THE LITHIUM-ION BATTERY BULGES SEVERELY, DISCONTINUE USE. KEEP SMALL CELLS WHICH ARE CONSIDERED SWALLOWABLE OUT OF REACH OF CHILDREN. SWALLOWING MAY LEAD TO BURNS, PERFORATION OF SOFT TISSUE, AND DEATH. SEVERE BURNS CAN OCCUR WITHIN 2 HOURS OF INGESTION. IN CASE OF INGESTION OF A CELL, SEEK MEDICAL ATTENTION PROMPTLY.

IMPORTANT SAFETY INSTRUCTIONS:
WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE TAKEN, INCLUDING THE FOLLOWING: READ ALL INSTRUCTIONS BEFORE USING APPLIANCE. SAVE THESE INSTRUCTIONS USE APPLIANCE ONLY FOR INTENDED USE AS DESCRIBED IN THIS INSTRUCTION BOOK. THE USE OF ATTACHMENTS NOT RECOMMENDED BY STYLPRO™ MAY CAUSE HAZARDS.

1. Please turn the device OFF after use. 2. This product should not be used by; children, pregnant women, patients with epilepsy, tumors and cancer patients, patients with cardiac pacemakers or similar. 3. Do not place any part of the product near heat sources such as radiators, cookers, etc do not lay on a soft surface such as a bed as device may overheat. 4. Keep the device away from corrosive liquids. 5. If you feel abnormalities in your skin during use or after use, please consult your medical specialist. 6. Do not tamper with any internal components, there are no user serviceable parts inside the unit. Do not attempt to repair the product yourself. 7. Keep away from water. This product is not a toy. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision. Please store with care and out of reach of children and pets. 8. To protect against electrical hazards, do not use while bathing or in the shower. 9. Close supervision is necessary when any appliance is used by or near children or persons made weak or disabled by illness or injury. 10. Do not use on sensitive or broken skin. Do not use on swollen or sunburnt skin or if you have any other skin conditions. 11. Remove all metals and solid objects, such as earrings, necklaces, bracelets, watches etc before use. Remove, glasses, false eyelashes and contact lenses before use. 12. Do not operate this appliance if it malfunctions or has been dropped on the floor or in water or damaged in any manner. 13. Do not directly touch the eyeball. 14. Do not share this product with anyone else. 15. Please do not use during charging.

LIMITED TWO-YEAR WARRANTY: OUR GOODS COME WITH GUARANTEES THAT CANNOT BE EXCLUDED UNDER THE EUROPEAN CONSUMER LAW. EU DIRECTIVE 1999/44/EC, THIS DIRECTIVE STATES A TWO-YEAR GUARANTEE APPLIES FOR THE SALE OF ALL CONSUMER GOODS EVERYWHERE IN THE EU. YOU ARE ENTITLED TO A REPLACEMENT OR REFUND FOR A MAJOR FAILURE YOU ARE ALSO ENTITLED TO HAVE THE GOODS REPAIRED OR REPLACED IF THE GOODS FAIL TO BE OF ACCEPTABLE QUALITY AND THE FAILURE DOES NOT AMOUNT TO A MAJOR FAILURE. STYLPRO™ WILL REPAIR OR REPLACE (AT OUR OPTION) YOUR APPLIANCE FREE OF CHARGE FOR TWENTY-FOUR MONTHS FROM THE DATE OF PURCHASE; IF THE APPLIANCE IS DEFECTIVE IN WORKMANSHIP OR MATERIALS. THIS WARRANTY IS AN ADDITIONAL BENEFIT AND DOES NOT AFFECT YOUR LEGAL RIGHTS. TO OBTAIN SERVICE UNDER THIS WARRANTY, PLEASE READ THE WARRANTY LEAFLET INCLUDED OR CONTACT STYLPRO™ CUSTOMER SERVICE AT TEAMTOM@STYLIDEAS.COM +44 (0) 20 3225 5559. TO MAKE A WARRANTY CLAIM YOU MUST RETAIN YOUR PROOF OF PURCHASE.

CORRECT DISPOSAL OF THIS PRODUCT: WHEN THE PRODUCT HAS REACHED THE END OF ITS USAGE, PLEASE DO NOT DISPOSE WITH REGULAR HOUSEHOLD WASTE. RECYCLING IS THE SAFEST AND MOST RESPONSIBLE WAY TO DISPOSE OF ELECTRICAL ITEMS THAT COME TO THE END OF THEIR LIFECYCLE. PLEASE DISPOSE OF IT ACCORDING TO YOUR LOCAL REGULATIONS. THE WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE) DIRECTIVE IN EUROPEAN COUNTRIES STIPULATES CORRECT RECYCLING PRACTICES. DO NOT DISPOSE IN REGULAR HOUSEHOLD WASTE.

CE MARKING: UK/EU -DECLARATION OF CONFORMITY: STYLIDEAS LTD DECLARES THAT STYLPRO COMPLIES WITH THE ESSENTIAL REQUIREMENTS AND OTHER RELEVANT PROVISIONS. A COPY OF THE DECLARATION OF CONFORMITY IS AVAILABLE ON REQUEST. STYLIDEAS LTD, AMSHOLD HOUSE, GOLDINGS HILL, LOUGHTON, ESSEX, IGIO 2RW, UK. OBELIS S.A, BD. GENERAL WAHIS 53, B-7030 BRUSSELS, BELGIUM.

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