
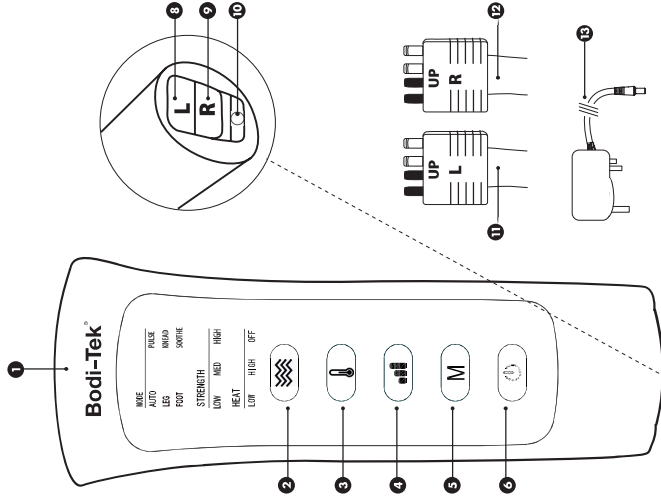
Comfort360° Air Compression Leg Massager

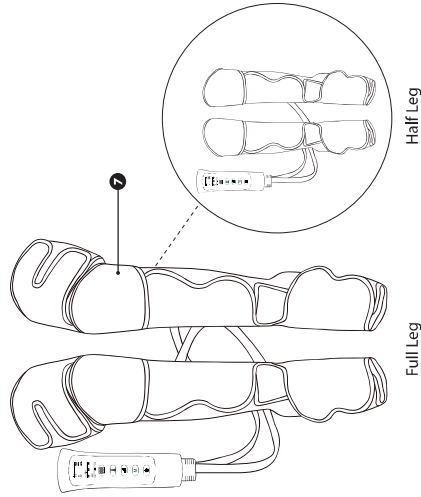Bodi-Tek®

Warnings

- Read all instructions thoroughly and retain for future reference.
- This appliance can be used by children aged 8 years and above and persons, with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance using adapter PSU (IP: 120V-1500B supplied).
- If the power lead or AC/DC adapter are damaged, it must be replaced by the manufacturer or its service agent or similarly qualified person in order to avoid a hazard. Contact Customer Care to arrange a replacement.
- Keep out of reach from children and pets.
- If you suffer from a skin condition, consult your dermatologist before use.
- Do not use on areas that have recently had surgery.
- Do not leave the device unattended. Always unplug when not in use.

Do NOT insert pins.

NOTE: These instructions cover both, the Comfort360 Leg & Full Leg Compression Massagers.

Know your unit

- 1 Controller
- 2 Vibration
- 3 Heat
- 4 Strength
- 5 Mode/Program (Auto, Leg, Foot, Pulse, Knead, Soothe)
- 6 On / Off
- 7 Hook & Loop Straps
- 8 Upper Socket
- 9 Lower Socket
- 10 Power Adapter Socket
- 11 Left Leg Connector
- 12 Right Leg Connector
- 13 Power Adapter

Getting Started

Ensure you are in a seated position. Use the Hook and Loop Straps to adjust the fit. Only a loose fit is required.

Plug both the Left and Right Lead Connectors into the Controller.

Insert the Right Leg Connector into the Lower Socket on the Controller with the UP facing upwards.

Insert the Left Leg Connector into the Upper Socket on the Controller with the UP facing upwards.

Push both in firmly.

Connect the Power Adapter to the Power Adapter Socket and plug into the mains supply.

How to use

- Press and hold the Power Button for 2 seconds to turn on. The boots will start to compress.
- Note: The Leg Massager has a memory function and will restart on the last used settings.
- Note: Press and hold the Power Button for 2 seconds to turn off.

Press the M / Mode/Program to select the following:

AUTO - Simultaneous compression pulsing through feet and legs.

LEG - Leg only compression.

FOOT - Foot only compression.

PULSE - A gradual squeezing sensation that slowly intensifies

KNEAD - A higher level of compression that slowly intensifies, starting at the feet and up to the legs

Press the Vibration symbol to change between:

Constant (White)

Intense (Green)

Intermittent (Blue)

You will feel the Vibration on the calves.

Select from: Low / High / Off

High Heat - 49°C

Low Heat - 38°C

The Leg Massager will automatically switch off after 20 minutes.

After Use

- The skin may appear reddened after use. This should fade quickly.
- Your muscles may also feel a little sore. This is normal and will lessen with regular use.
- Tip: Ensure you stay hydrated to help flush out any toxins and improve muscle recovery time.

Cleaning & Maintenance

- Switch off and clean using a soft damp cloth.
- Do not immerse the device/unit in water or other liquids.
- Dry thoroughly before storing.

Important Information

There are no user serviceable parts inside the appliance.

Disposal

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic Equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:


1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly

Specifications

Model: BT-MAFL/BT-MAHL
Input: 12V 1500mAh 18W

Read these instructions carefully before use and then retain for future use.

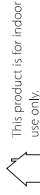This product is for indoor use only.

This unit complies with current EU product legislation.

Centre pin positive

Extension Straps

Extension Straps for a larger fit are available at:

www.bodi-tek.com/collections/spare-parts

Guarantee

To register your two year guarantee online, visit:

www.bodi-tek.co.uk/warranty

Customer Care

If you have any questions, please contact us at: online@bodi-tek.co.uk

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© The Dezac Group Ltd 2023, PO Box 17, Cheltenham Spa, England GL53 7ET.
© The Dezac Group Ltd 2023, PO Box 17, Cheltenham Spa, Ingilterra GL53 7ET.
RP Regno Unito: The Dezac Group Ltd, PO Box 17, Cheltenham Spa,
EU RP: The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB