

WATER FLOSSER & ORAL IRRIGATOR

the beauty specialists

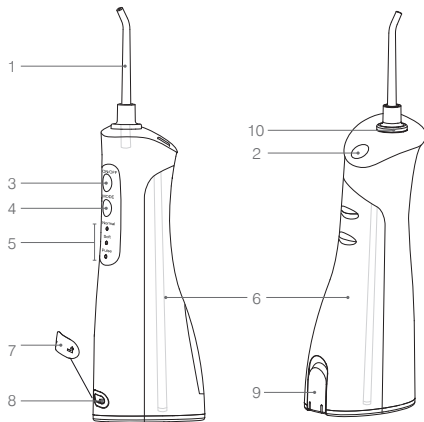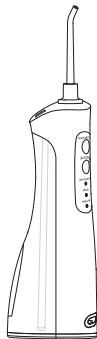ENGLISH	4
FRANÇAIS	14
DEUTSCH	24
NEDERLANDS	34
ESPAÑOL	44
PORTUGUÊS	54
ITALIANO	64

WARNINGS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance. If any part of the mains adapter or cord are damaged, discontinue use. Contact Customer Care to arrange a replacement.
- Do not immerse product in water. Do NOT use solvents or abrasive cleaners.

- For external use only.
- Do NOT use this product if it is not working properly, or if it has been dropped or damaged.
- Ensure the unit is off and disconnected from the power source when filling with water and cleaning the water reservoir.
- Keep out of reach from children and pets.
- Do not store the product where it can fall or be pulled into the bath or sink.
- Keep away from any appliances which are sensitive to water to avoid accidents.
- To prevent bacteria from propagating, ensure the unit and nozzle(s) are thoroughly washed with soapy water and leave to dry before storing.
- Do not fill the reservoir with water warmer than 40°C.

- Fill the water reservoir with water or other dental professional recommended solution only (i.e. mouthwash or other prescribed solutions). Some antibacterial solutions can shorten the life of the product.
- Do not use for any other purpose than oral cleaning.
- If you have had oral or dental surgery in the past 2 months, consult your dentist before using this product.
- Do not press the nozzle too firmly against your teeth or gums. Doing so may cause injury.
- We recommend replacing your nozzle every 6 months.
- Do not use if you have an open wound on your tongue or in your mouth.
- Do not use while wearing any oral jewellery. Remove prior to use.
- Do not use with mouthwashes which contain iodine, bleach or tea tree oil.

YOUR UNIT

- | | |
|--------------------------------------|------------------------------|
| 1 Nozzle | 6 Water reservoir |
| 2 Nozzle release button | 7 Charging socket cover |
| 3 ON/OFF button | 8 Charging socket |
| 4 Programme (MODE) button | 9 Water reservoir supply cap |
| 5 Programme (MODE) indication lights | 10 Nozzle socket |

IMPORTANT: There may be some liquid remaining inside the water reservoir prior to first use. This is the remaining distilled water which was used for the water flow performance test, so there are no health or hygiene concerns.

CHARGING

IMPORTANT: Fully charge the unit for 24 hours before using for the first time. Charging will take approximately 4 hours thereafter. Plug the USB cable into the charging socket on the unit and connect to a power source.


Important: Ensure there is no water remaining inside the socket before connecting to the power source.


The 'Pulse' mode indicator light will illuminate white whilst the unit is charging. It will illuminate white once the unit is fully charged.

SET UP

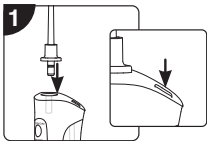Insert the nozzle into the nozzle socket.

Press the nozzle release button to remove the nozzle.

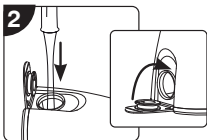Open the water reservoir supply cap and fill with **lukewarm** water.

Click shut.

Important: Ensure the unit is turned OFF when filling up the water reservoir.

PROGRAMMES (MODE)


Normal – for removing particles of food trapped between teeth.

Soft – for sensitive gums.

Pulse – for deeper penetration, removing plaque causing bacteria, massaging gums and rinsing between teeth and gums.


Lean over a sink and keep your mouth **slightly** open to allow the water to circulate and flow out. Due to water spray during treatment we advise using over a sink.

Note: You may need to refill the water reservoir during treatment.

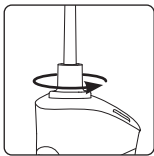You can direct the water stream to target hard to reach areas by turning the 360° nozzle attachment.


USE

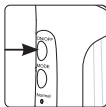
Aim at the teeth and press the ON/OFF button to turn the unit ON. Press the Programme (MODE) button to select a programme and start irrigating.

Note – the unit will automatically turn off after 2 minutes.

Direct the water jet flow at a 90° angle to the tooth at the gumline and work slowly along the teeth and gum line cleaning inside and outside of both the upper and lower teeth.

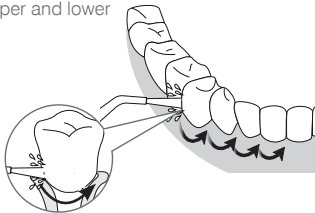Clean thoroughly between teeth and around any braces, crowns or bridgework.

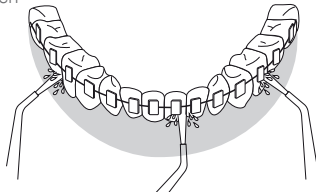Stop using the product if you experience pain or sensitivity during treatment and seek professional advice from your dentist.

AFTER EACH USE

Turn the unit OFF.

Open the water reservoir supply cap and empty any remaining water away.

Press the ON/OFF button to drain any liquid still remaining inside the unit and nozzle.

Press the nozzle release button to remove the nozzle.

Important: To prevent bacteria from propagating, ensure the unit and nozzle(s) are thoroughly washed with soapy water and leave to dry before storing.

SPECIFICATION

USB CHARGER RATED INPUT

DC 5V /500mA

READ THESE INSTRUCTIONS CAREFULLY BEFORE USE AND THEN RETAIN FOR FUTURE USE.

THIS PRODUCT IS FOR INDOOR USE ONLY.

THIS PRODUCT COMPLIES WITH CURRENT EU PRODUCT LEGISLATION.

Model: DCIR2

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly

CUSTOMER CARE

If you have any questions, or if you experience any difficulties please call our Careline on: **01242 702345** or email us at: **online@riobeauty.com** To reorder additional accessories, please visit **www.riobeauty.com** or contact your local distributor.

GUARANTEE

To register your two year guarantee online, visit: **www.riobeauty.com/warranty_online.htm**

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© The Dezac Group Ltd 2022, PO Box 17, Cheltenham Spa, England GL53 7ET.
RP REGNO UNITO: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, Inghilterra GL53 7ET
RP UE: Il gruppo Dezac (Europa) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Paesi Bassi
© The Dezac Group Ltd 2022, PO Box 17, Cheltenham Spa, Inghilterra GL53 7ET.
UK RP: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, England GL53 7ET
EU RP: The Dezac Group (Europe) B.V. Azura,