AILORIA

PRO SMILE

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INSTRUCTION MANUAL

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Schallzahnbürste USB

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Sonic Toothbrush USB

Product features

- Up to 40,000 brush head movements per minute
- 5 cleaning modes:
Sensitive - Clean - Whitening - Massage - Polish
- 2 minute timer and 30 second interval for 4 quadrants
- Charging of the lithium-ion battery via USB
- Charging time approx. 4 hours
- 25 days battery life with 2x daily use of 2min
- Waterproof according to IPX6

Getting started

Charge the toothbrush before using it for the first time.

Charging time is approximately 4 hours to charge the toothbrush before first use or fully discharged battery.

Put the brush head on the toothbrush and make

sure it is locked in place.

Charging

If the red indicator light flashes or the product does not turn on, the battery power is insufficient and the product needs to be charged.

To charge, plug the supplied cable into the bottom of the toothbrush and connect the other side to a USB port or charging adapter.

To prevent overcharging the battery, charging will automatically stop when fully charged and the indicator light will turn off.

When using for the first time or if the battery has been discharged for a longer period of time, the indicator lights up for the first 5 minutes.



Installing the brush head

Put the brush head on the neck of the toothbrush. Align the guide mark on the back of the brush head with the beginning of the guide mark on the back of the brush body.

Rotate the brush head so that the guide marks on the back of the brush head and neck are aligned.

Dentists recommend replacing the brush heads

every 3 months to achieve optimum results.

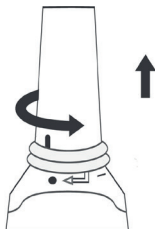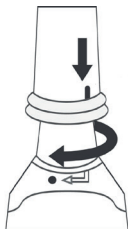Application

1. Switch on the toothbrush.
2. During the entire brushing time, move the brush head slightly inclined (approx. 45°) and slowly without pressure over the teeth and chewing ridges. The brush head should be placed with one half on the gums and the other half on the teeth to ensure cleaning even on the important transition between gums and teeth.
3. Brush your teeth for 30 seconds on each of the quadrants. The Interval Timer indicates when the quadrant should be changed by briefly interrupting. The Smart Timer ensures the two-minute brushing time recommended by the dentist.

NOTE: PRO SMILE is safe to use on braces and dentures (fillings, veneers, crowns, bridges, etc.). When used on braces, the brush heads wear out more quickly.

Dentists recommend replacing brush heads every 3 months for best results.

PRO SMILE has been tested and meets the safety standards for electromagnetic devices.

Cleaning modes

The PRO SMILE sonic toothbrush has 5 cleaning modes: Sensitive - Clean - Whitening - Massage - Polish. The memory function remembers which mode you used last.

Press the power button and switch on the toothbrush. Press the power button again within 3 seconds to switch to the next mode. Repeated pressing activates all 5 modes consecutively. After the fifth and last mode, the toothbrush switches off.

The indicator lamps on the toothbrush show you which mode is activated. When you switch to the first mode „Sensitive“, the first indicator lamp lights up. When you switch to the 2nd mode „Clean“, the second indicator lamp lights up and in the 3rd mode „Whitening“, the third indicator lamp lights up. When you switch to the 4th mode „Massage“, the first and

second indicator lamps light up simultaneously. In the 5th and last mode „Polish“, the first, second and third indicator lamps light up simultaneously.

Mode	Display	Application
Cleaning modes		
Sensitive		1. indicator lamp lights up For cleaning sensitive teeth and gums
Clean		2. indicator lamp lights up For daily cleaning of normal teeth and gums
Whitening		3. indicator lamp lights up For particularly intensive cleaning of the tooth surfaces with whitening effect
Additional modes		
Massage		1. and 2. indicator lamp lights up For gum massage with gentle vibrations of the brush head
Polish		1., 2. and 3. indicator lamp lights up For extra treatment to remove discolouration and polish tooth surfaces

Cleaning

Do not clean the PRO SMILE sonic toothbrush and brush heads in the dishwasher.

Rinse the brush head and bristles with water after each use.

Remove the brush head and rinse the neck of the toothbrush with warm water.

Use a clean, damp cloth to wipe the entire surface of the toothbrush.

If you do not use PRO SMILE for a long period of time, remove the charging cable, clean the toothbrush and store it in a cool, dry place and avoid direct sunlight. It is recommended to recharge the toothbrush after 6 months of non-use.

Warnings

- Consult your doctor before using PRO SMILE if you have any medical problems.
- If you have had jaw or gum surgery in the last two months, consult your doctor before using PRO SMILE.
- If excessive bleeding occurs after using PRO SMILE or if bleeding continues after one week of use, consult your dentist.
- Do not clean PRO SMILE or the brush heads in the dishwasher.
- Do not use attachments other than those recommended by the manufacturer.
- Discontinue use if PRO SMILE appears damaged in any way (brush head, handle).
- Keep the cable away from all sources of heat.

- Do not use a brush head with damaged bristles. Deformed bristles can break off when brushing. Replace the brush head every 3 months or sooner if signs of wear appear. Be aware that improper cleaning can damage the bristles. (See Cleaning)
- Use PRO SMILE only for its intended use as described in this manual. This product is designed to clean teeth and gums. If you experience pain or discomfort, please stop using PRO SMILE and consult your dentist or doctor as this may indicate other health problems.
- Persons (including children) with lack of experience and knowledge or reduced mental, sensory or physical capabilities should not use this product unless they are supervised during use or unless they have been given appropriate instructions concerning the use of PRO SMILE by a person responsible for their safety.
- If your toothpaste contains peroxide, baking soda or other bicarbonate, clean your brush head and treat it thoroughly with water after each use. This will prevent plastic breakage due to prolonged exposure to these chemicals.

Troubleshooting

The toothbrush will not switch on

You may need to charge PRO SMILE.

Tingling sensation when brushing

Due to the sonic vibrations, you may experience a tingling sensation when using PRO SMILE for the first times. This feeling will decrease over time as you get used to the vibrations.

The sonic vibrations seem weak

The bristles may be pressed too hard against your teeth. Gently place the brush head on your teeth.

Warranty & Service

Defects due to material and/or workmanship will be repaired or replaced by AILORIA under warranty if proof of purchase is provided during the warranty period. Contact the Customer Service Centre at service@vehnsgroup.com.

Items not covered by the warranty:

- Brush heads
- Damage caused by misuse, abuse, neglect, alteration or unauthorised repair
- Normal wear and tear including scratches, discolouration or fading

Disposal

According to the EU directive, this product must not be disposed of together with other household waste. Find out about and follow local regulations on the separate collection of electrical and electronic equipment. Proper disposal of old products helps to prevent negative effects on the environment and health.

Manufacturer

VEHNS GROUP GmbH
Theatinerstrasse 40-42
80333 Munich
Germany

www.vehnsgroup.com
service@vehnsgroup.com
www.ailoria.de

MADE IN CHINA

Changes and errors excepted

The manufacturer accepts no liability for damage caused by incorrect, improper or incompatible use of the product.

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© VEHNS GROUP GmbH
VEHNS GROUP GmbH, Theatinerstrasse 40-42