
SHINE BRIGHT

DE Bedienungsanleitung **EN** User manual **FR** Manuel de l'utilisateur **IT** Manuale di utilizzo **ES** Manual de Usuario **NL** Gebruiksaanwijzing **PL** Instrukcja obsługi **LT** Naudotojo vadovas **CS** Návod k použití **DA** Instruktionsmanual **FI** Käyttöohje **NO** Bruksanvisning **SK** Používateľská príručka **SV** Bruksanvisning **HR** Korisnički priručnik **PT** Manual do utilizador **BG** Ръководство за потребителя **RO** Manual de utilizare

Version 1.0 - 05/2025

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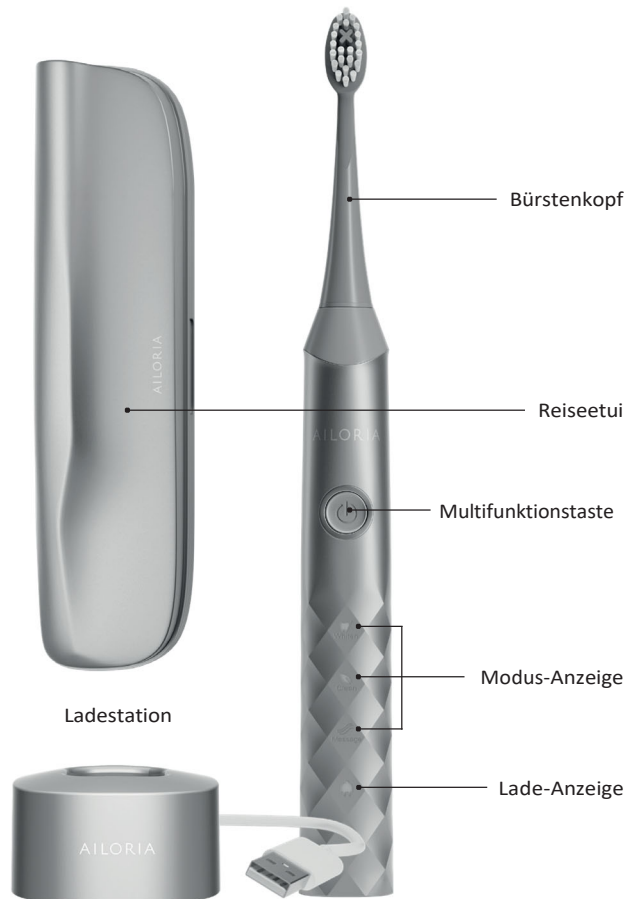
Sonic Tooth Brush

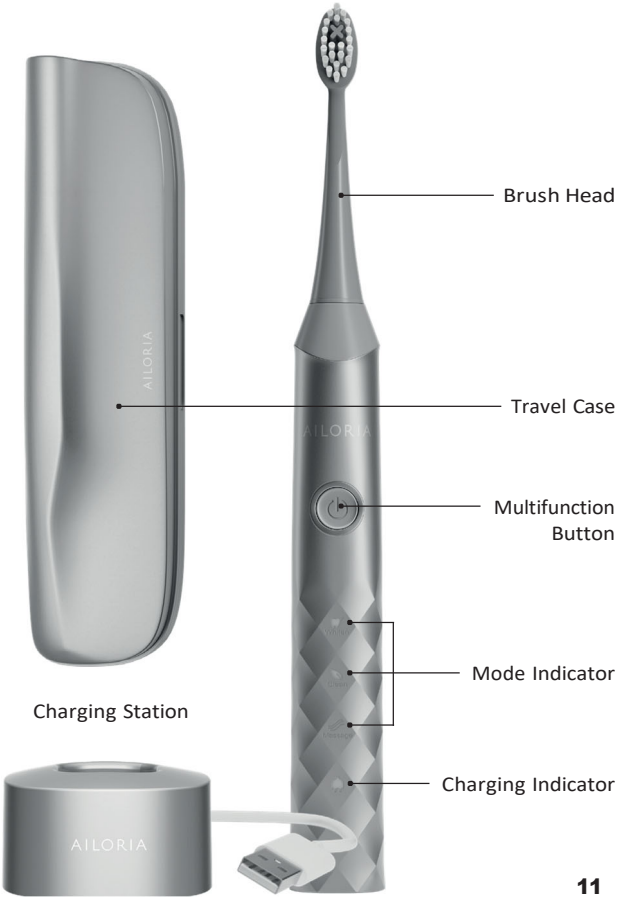PRODUCT FEATURES

- Modes: Whitening, Clean, Massage
- Timing: 2-minute timer with 30-second intervals
- Battery: 700 mAh lithium-ion, up to 21-day life
- Charger: Inductive USB, 1.5W/DC 5V, with safety features
- Durability: IPX7 waterproof standard

SAFETY INSTRUCTIONS

- Use the sonic toothbrush only as directed in this manual.
- If you've had dental surgery in the past two months, consult with your dentist before using.
- Stop use and consult your dentist if you experience excessive gum bleeding, discomfort, or pain.
- This toothbrush should not interfere with pacemakers, but consulting a healthcare provider is recommended for safety.
- Children under 8 and individuals with reduced physical, sensory, or mental capabilities should use this toothbrush only under supervision.
- Keep the toothbrush out of reach of children at all times.
- Do not use the toothbrush if the brush head, handle, charger, or USB cable is damaged.
- Repairs should be carried out by the manufacturer only.
- Avoid using damaged brush heads; replace every 3 months or sooner if signs of wear are visible.
- Use only manufacturer-recommended attachments. Do not clean the toothbrush, attachments, or charger in the dishwasher.
- Clean the toothbrush under running water; improper cleaning may damage the bristles.
- Clean the charger with a damp cloth only; it should not be

washed under running water. Unplug before cleaning.

- If using toothpaste with peroxide, baking soda, or bicarbonate, rinse the brush head thoroughly after each use to prevent damage.
- Keep the charger and USB cable dry and away from liquids.
- Ensure the device's voltage matches your local voltage before connecting.
- Check the compatibility of the mains plug with the product's USB cable.
- The USB cable is permanently attached and should be replaced only with an original part if damaged. Contact your dealer for assistance.
- Avoid excessive bending of the USB cable to preserve its function.
- Do not place the charger on hot surfaces.
- Avoid overcharging the battery; disconnect the device from the charger once it is fully charged.
- Store the device in a cool, dry place and steer clear of extreme temperatures.
- Do not expose the battery to direct sunlight and do not store the device in damp environments.
- Protect the device from impacts, drops, and excessive pressure to prevent damage to the battery.
- Do not attempt to remove or replace the battery yourself.
- If you notice any signs of damage or leakage from the battery, stop using the device immediately and contact customer service.
- Avoid contact with leaking fluids; if skin contact occurs, wash thoroughly with water right away.
- Dispose of the device at the end of its life cycle according to local regulations for electronic waste.

- Do not throw the battery into fire or water.
- If the device overheats, emits smoke, or shows any other signs of malfunction, turn it off immediately and disconnect it from any power source.
- Keep the device away from flammable materials until the issue is resolved.

CHARGING

Before first use, charge the toothbrush for approximately 8 hours, or until fully charged. A flashing red charging indicator or inability to power on signals a need to recharge.

If the charging plug symbol flashes red or the product cannot be switched on, the battery power is insufficient and the product must be recharged. Connect the charger to a power source using the USB adapter provided, or any standard 5V/1A USB power supply if not included. Charging can also be done via other USB sources, like a laptop or PC. The charging begins automatically once placed on the charging station, with the symbol lighting up red. The indicator turns off when fully charged, offering up to 21 days of use (based on twice-daily, 2-minute sessions). The toothbrush is inoperable while charging and will automatically switch off and enter charging mode if placed on the station while on.

After charging, disconnect the charger and store it safely. Clean with a damp cloth only; do not submerge in water.

INSTALLING THE BRUSH HEAD

Slide the brush head onto the metal shaft until it clicks into place securely. To remove it, simply pull the brush head straight up. Avoid using brush heads with damaged or deformed bristles, as they can break off during use. It's recommended to replace the brush head every 3 months.

CLEANING MODE

Initially, powering on the sonic toothbrush activates the Whitening mode. Quickly press the On/Off button again within 3 seconds to cycle to the next mode. Continue pressing to cycle

through all three modes: Whitening, Clean, and Massage. The toothbrush turns off after the final mode. Indicator lights display the active mode. The toothbrush features three modes:

Whitening: Intense cleaning and plaque removal, best for those without sensitive teeth.

Clean: Gentle daily cleaning, suitable for sensitive teeth.


Massage: Varying intensity to stimulate gums and polish teeth.

To turn off the toothbrush, press the On/Off button 3 seconds after the last activation. Pressing too soon cycles modes. The toothbrush remembers the last mode used and will start in that mode next time.

HOW TO USE

Wet the brush head, apply toothpaste, and turn on the toothbrush by pressing the on/off button once it's in your mouth.

Throughout brushing, angle the brush head about 45° to your teeth and gum line, moving it slowly and gently across your teeth and molars without applying pressure. Position the brush so that half is on your gums and the other half on your teeth, targeting the crucial junction between gums and teeth for optimal cleaning. The premium bristles, enhanced by the sonic technology, ensure a thorough yet gentle clean.

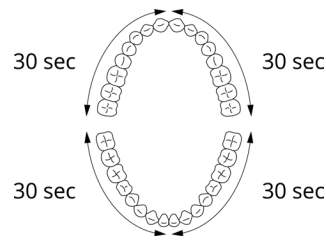Brush each quadrant of your mouth for 30 seconds. The interval timer pauses briefly to signal when to move to the next quadrant. A smart timer guarantees you meet the dentist-recommended two-minute brushing time.

CLEANING & STORAGE

After each use, rinse the brush head and bristles under running water. The SHINE BRIGHT toothbrush is waterproof due to its IPX7 rating, allowing for shower use—just keep the charger dry. Wipe the handle with a damp cloth.

Weekly, detach the brush head to rinse the metal shaft with warm water. Avoid damaging the rubber seal with sharp objects. For best performance, replace brush heads every 3 months or upon signs of wear, using only the manufacturer's recommended replacements.

Before cleaning the charging station, unplug it. Wipe the charger and its cover with a damp cloth. Do not submerge the charger or expose it to running water.

For long-term storage, disconnect the USB cable or plug, clean the toothbrush, and store it in a cool, dry place away from sunlight and children's reach. If storing the toothbrush for a longer period of time, recharge it every 6 months.

TROUBLESHOOTING

Toothbrush Won't Turn On: It might require charging.

Tingling Sensation: Initial use may cause a tingling due to the sonic vibrations. This sensation typically diminishes as you become accustomed to the vibrations.

Weak Vibrations: Applying too much pressure can weaken the perceived vibrations. Ensure you're moving the brush head gently over your teeth and gums.

ACCESSORIES & REPLACEMENT BRUSHES

Replacement brush heads, charging station and charging plug are available from your authorized dealer or directly from the manufacturer: www.ailoria.de

WARRANTY & SERVICE

This product comes with a two-year warranty starting from the date of purchase. The warranty does not cover:

- Brush heads

- Damage resulting from improper use, misuse, neglect, modification, or unauthorized repair
- Normal wear and tear, including scratches, discoloration, or fading

TECHNICAL SPECIFICATIONS

Toothbrush Power: 3 W
Toothbrush Voltage: DC 3.7 V
Battery Capacity: 700 mAh
Charger Power: 1.5 W
Charger Voltage: DC 5 V

DISPOSAL INSTRUCTIONS


Electrical and electronic equipment must not be discarded with regular household waste, but must be recycled. According to the Electrical and Electronic Equipment Act, such items are marked with a crossed-out trash bin symbol to encourage proper waste separation. You can return old electrical devices at designated free collection points, the addresses of which can be obtained from your city or local authority.

Batteries and accumulators that can be removed from the device without damage must be taken out before disposal. As a consumer, you are required to dispose of used batteries and accumulators at dedicated collection centers, ensuring they are completely discharged before disposal.

This product includes a rechargeable lithium ion battery (3.7V, 700mAh).

Subject to changes and errors

Specifications are subject to change without notice, and the manufacturer is not liable for damages resulting from incorrect, improper, or incompatible use of the product.

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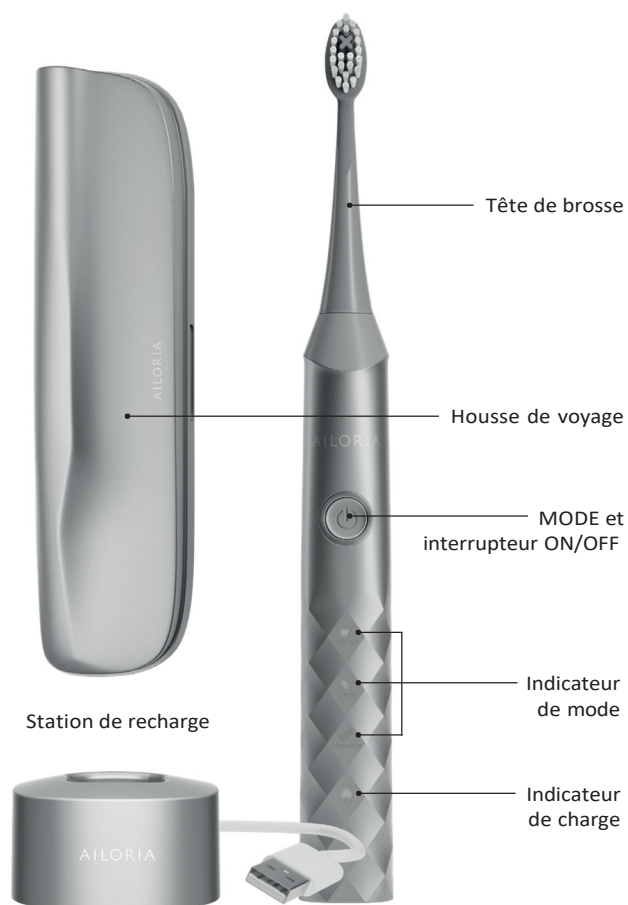A diagram of a circular exercise routine. It consists of a large circle divided into four equal quadrants by a vertical and a horizontal line. Each quadrant contains a different exercise icon: a person in a lunge position, a person in a squat position, a person in a lunge position, and a person in a squat position. Arrows on the outer edge of the circle indicate a clockwise flow from one quadrant to the next. Each quadrant is labeled '30 sec'.

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© VEHNS GROUP GmbH
HUNAN HUAHUI NEW ENGERGY CO.,LTD
VEHNS GROUP GmbH Theatinerstr. 40-42, 80333 Munich, Germany,
VEHNS GROUP GmbH