
SMART CONNECT

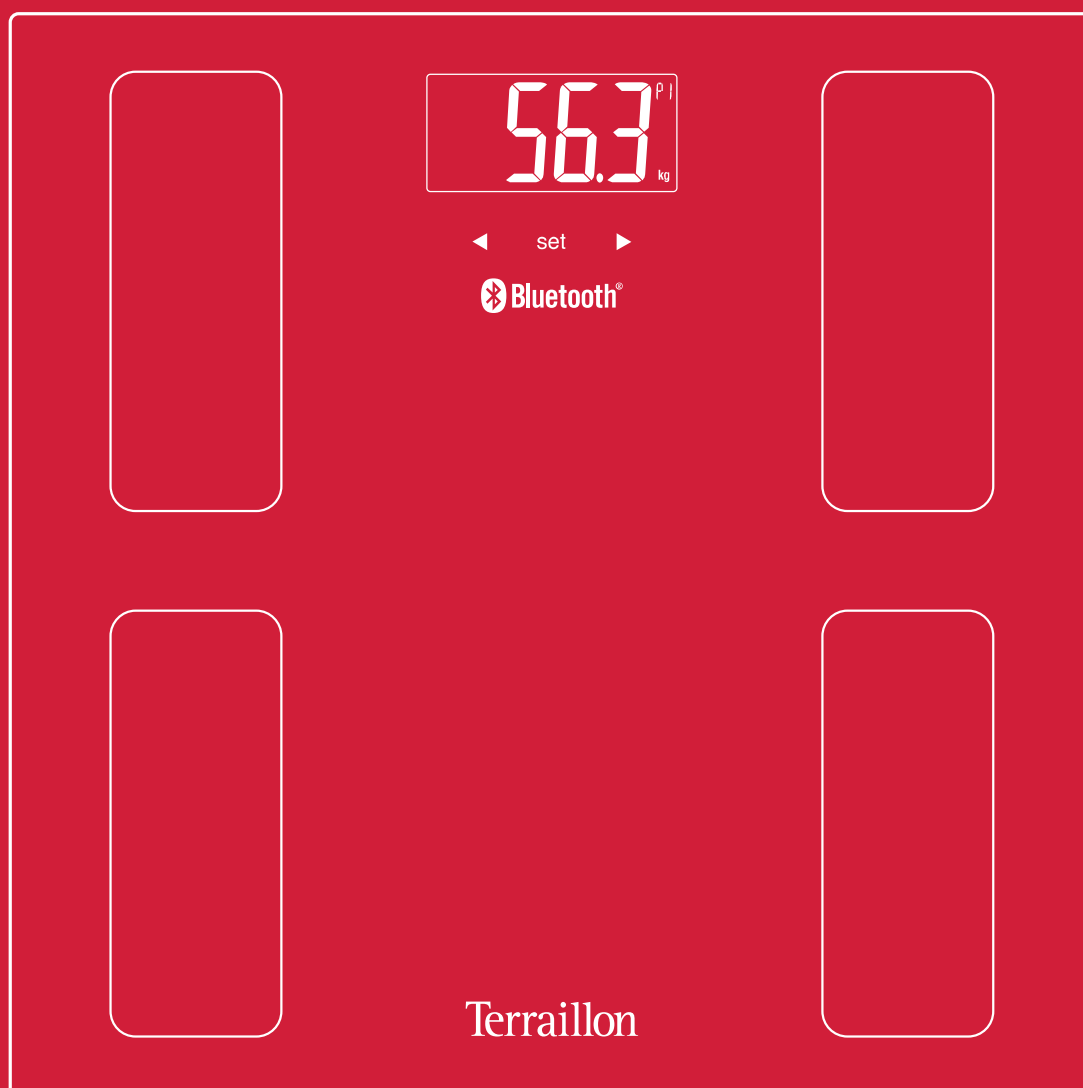USER MANUAL | GUIDE D'UTILISATION

NL Handleiding | **IT** Manuale d'uso | **ES** Manual del usuario

DE Bedienungsanleitung | **PT** Manual do utilizado

Terraillon®

www.terraillon.com

Made for / Compatibilité

iOS 8.0 & +

Android 4.3 & +
Smartphone only / No tablet

iPhone is registered trademark of Apple Inc.

Android is a trademark of Google Inc.

The Bluetooth® word mark and logos are registered trademarks
owned by Bluetooth SIG, Inc.

Other trademarks and trade names are those of their respective owners.

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Thank you very much for selecting the Terraillon Glass Body Fat Analyzer SMART CONNECT.



Don't use this product if you wear a pacemaker or implantable electronic device • Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening.

Need support?

<http://terraillonhelp.zendesk.com>

PRODUCT OVERVIEW

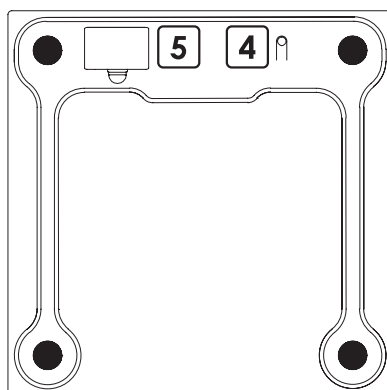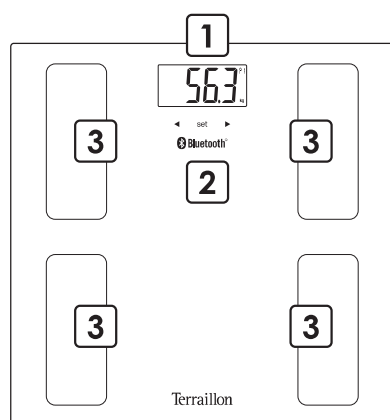- 1** LCD Display
- 2** Profile selection
- 3** Electrodes
- 4** Unit and pairing button
- 5** Battery compartment

FIRST USE

- Insert the batteries (3xAAA).
- You will need your smartphone to install your bathroom scales.
- Download the Wellness Coach – Health app, then open it.

- Follow the instructions that appear on your smartphone. Ensure you enter all the information requested.
- Your smartphone must remain near to the bathroom scales throughout the installation process.
- If you encounter any problems, please consult the help centre:

<http://terraillonhelp.zendesk.com>

- It's important to bind your account with the scale for automatic recognition of user during weighing. You can pair the scale with 8 smartphones maximum. To do that, just repeat steps B to F on an other smartphone and select another profile number. To bind an other user on the same smartphone, please log in with the other user account and repeat steps E and F. Then select another profile number.

Before you first weighing, please select your profile number on the scale.

When the scale does not recognize user, two profiles are displayed (here P1 and P8). Press ◀ to select profile on the left, ▶ for the one on the right.

START MEASURING

1. You must have bare feet to be able to analyse your body mass.
2. Open the "Wellness Coach Health" application.
3. Step onto your bathroom scale. Maintain contact with the electrodes until the following progress bar disappears from the screen.

4. If the application is not open, the data will be saved on the bathroom scale.
5. The symbol  will then be displayed on the screen.
6. To send the data, open the application on your smartphone, select your profile on the scale, then press SET. Otherwise, it will be sent the next time you weigh yourself if you follow step A.

Be sure to always weigh yourself on a hard and flat surface.

MAINTENANCE

Use water only. DON'T use propellant or other chemicals.

WARRANTY & ENVIRONMENT PROTECTION

This product is guaranteed for 2 years against material and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase must be presented if the scale is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint please contact the store where you purchased the scale first.

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

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Terraillon SAS France & Headquarters
Terraillon UK Limited
Terraillon Asia Pacific Ltd
Terraillon Corp USA