

RADIANCE LIGHT THERAPY PANEL

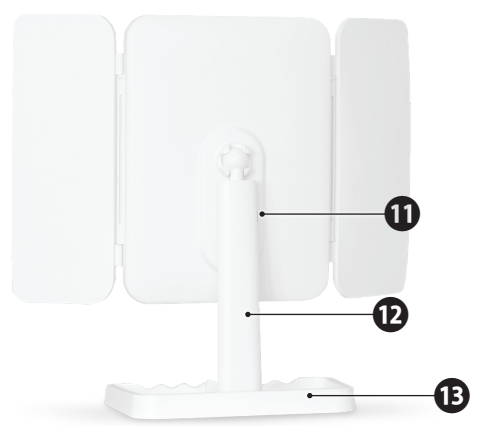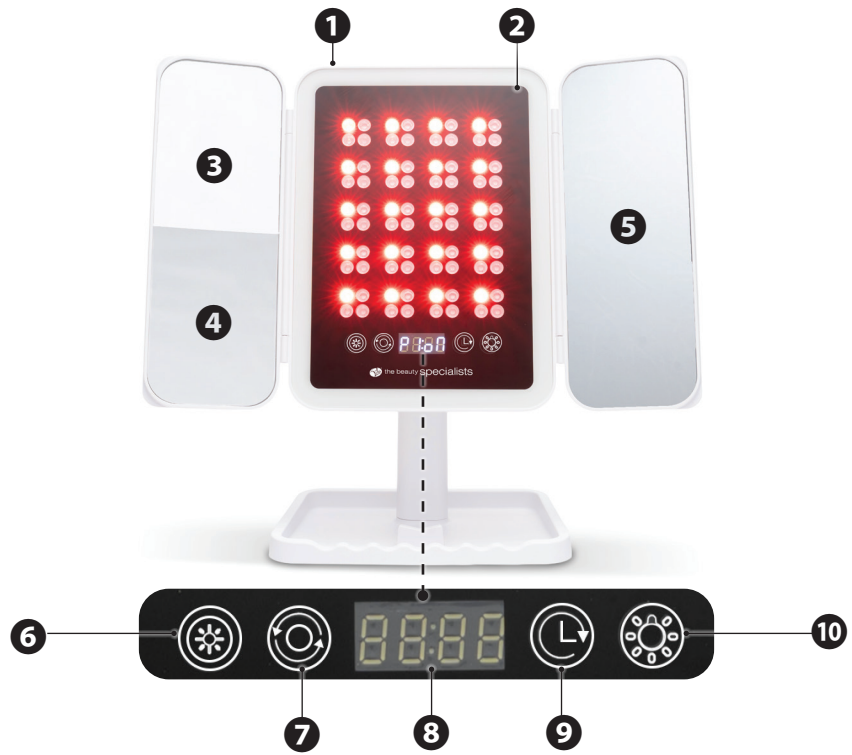CONTENTS


ENGLISH	4
FRANÇAIS	10
DEUTSCH	16
NEDERLANDS	22
ESPAÑOL	28
PORTUGUÊS	34
ITALIANO	40

WARNINGS & PRECAUTIONS

- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Keep away from direct sunlight when not in use. The magnification mirrors could cause a fire.
- **Please read all instructions thoroughly and retain for future reference.**
- Keep out of reach from children and pets.
- This appliance can be used by children aged 18 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Do not place or store in direct sunlight due to the risk of fire from reflected/magnified light. Store in a cool dry place.
- Do NOT use this appliance if it is not working properly, or if it has been dropped or damaged.
- Do not charge for longer than 24 hours.
- The USB-A to USB-C cable can be changed if damaged using a suitably approved USB cable.
- Do not leave the device unattended. Turn off before you leave the room.
- Do NOT immerse in water.
- The appliance may be cleaned by wiping with a damp cloth. Do NOT use solvents or abrasive cleaners.
- This device must only be used for its intended purpose as described in this user guide.
- Do not use the Light Therapy Panel if you are taking any medication that may cause sensitivity to light.
- Using the Light Therapy Panel while on light-sensitizing medication may lead to skin irritation or other reactions.
- Some medications, such as certain antibiotics, cancer treatments, and diuretics, can make your skin more sensitive to light, so always check with your doctor if you are unsure.
- Use the device only for the specific skin conditions listed in the instructions and not for any other purpose.
- Do not use the Light Therapy Panel if you are pregnant, breastfeeding, or trying to become pregnant.
- Avoid using the device if you have a medical condition like Lupus, Albinism, or photosensitive eczema.
- Do not use the Light Therapy Panel if you have any disorder that causes sensitivity to light.
- Always perform a patch test before using the device for the first time to check for any skin sensitivity.
- DO NOT use Light Therapy Panel in any way that is not exactly as specified in these instructions for use.
- Stop using the Light Therapy Panel immediately if it feels hot or causes your skin to feel uncomfortably warm.
- Make sure to clean your face and remove any makeup, creams, or lotions before using the Light Therapy Panel.
- Keep the device away from water and do not use it in wet environments or if it may have come into contact with liquids.
- Do not place the device near hot surfaces or expose it to high temperatures.
- Never try to modify or repair the device yourself as doing so can void the warranty and may be unsafe.
- Do not use the Light Therapy Panel if any part such as the Light Therapy Panel, or charging cable is damaged.
- Do not use the device if you are prone to headaches triggered by bright or flashing lights.
- Follow the recommended treatment guidelines and avoid using the Light Therapy Panel more often than advised.
- Do not use on areas that have had cosmetic surgery.
- If you suffer from a skin condition, consult your dermatologist before use; Sensitivity to light, or photosensitivity, is a known side effect of many medications. This may include certain types of antibiotics, cancer treatments, and water pills (diuretics). If you're not sure whether a medication you're taking could cause this reaction, it's best to check with your doctor or pharmacist.
- In addition to medications, other substances can also make the skin more reactive to sunlight. Some common examples include: St. John's Wort, coal tar, some deodorants, antibacterial soaps, artificial sweeteners, mothballs (which contain naphthalene), petroleum-based products, optical brighteners in laundry detergents, and cadmium sulfide (sometimes used in tattoo pigments).

INTRODUCTION

With aging, the skin gradually loses its firmness, moisture, and elasticity. Specific wavelengths of light have been shown to support skin rejuvenation and healthier-looking skin over time. The Radiance Light Therapy Panel uses this technology to offer a non-invasive treatment to help improve your skin's appearance.

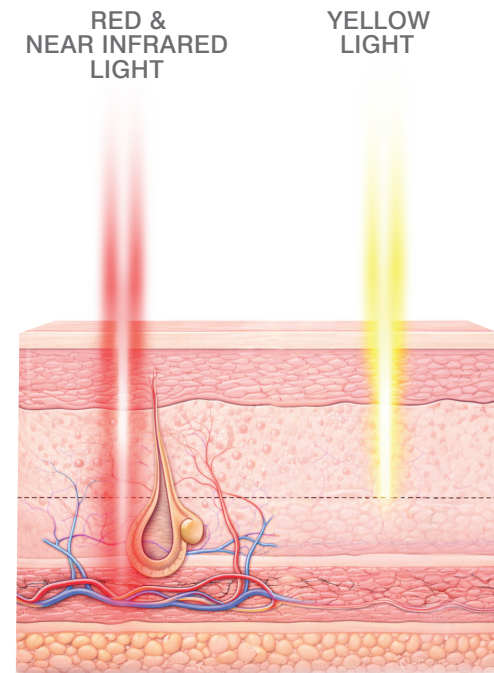
KNOW YOUR UNIT

- 1 Light Therapy Panel
- 2 Light (Cool, Natural, Warm)
- 3 2x Magnification
- 4 3x Magnification
- 5 1x Magnification
- 6 LED Therapy - On/Off & Intensity
- 7 Treatment Programme
- 8 Display (Timer, Mode & Charging)
- 9 Timer
- 10 On/Off/Light Colour (Cool, Natural, Warm)
- 11 USB Socket
- 12 Stand
- 13 Base
- 14 Charging Cable
- 15 Comfort Glasses

ASSEMBLY

1. Connect the Stand to the Light Therapy Panel aligning the three ribs and pushing firmly into place.

2. Push the Stand into the Base. The Stand will only push into place when the small locating points are aligned.

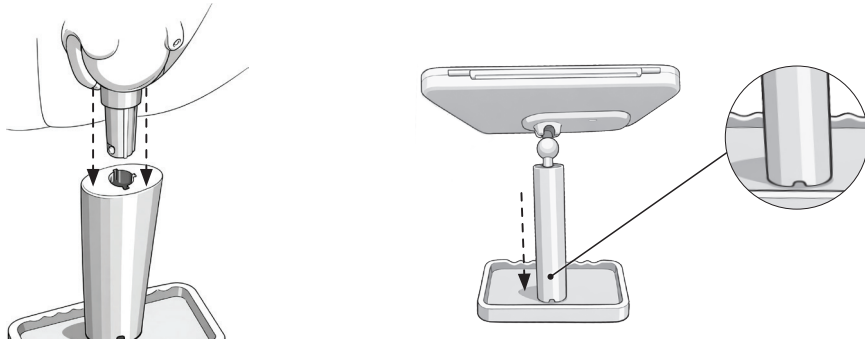
PATCH TEST

Before your first treatment, it's important to conduct a sensitivity test: Place unit over your inner forearm and run the Treatment Programme and intensity level you intend to use for three minutes. If you experience redness or itching which persists for more than an hour do not continue to use the product.

CHARGING

Always recharge when the display shows "LO" for low battery.

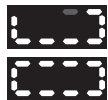To charge, insert the USB Charging Cable and plug into a USB port. It will take between 3-3.5 hours to fully charge.

PREPARING FOR USE

1. Remove makeup and other skincare products.
2. Put on the Comfort Glasses.

GETTING STARTED

It is recommended to use the Light Therapy Panel consistently over several weeks. Start with daily sessions for 3-4 weeks. It is optimal to have 15 minutes of treatment per day. This can be done in a number of shorter intervals, for example, three 5 minute sessions. After the initial 3-4 weeks the frequency may be reduced to 3 to 5 times a week, or to suit your skin condition or personal skincare objectives.

TREATMENT OPTIONS

Treatment	Distance from LED Panel	Treatment Time	Use
Intense/Targeted	5-10 cm	5-15 minutes	Focused treatment while sitting upright
Relaxed/Comfort	10-30 cm	15-60 minutes	Gentle exposure while relaxing or reading

 The LEDs emit high intensity light. For your comfort the glasses should be worn as you bring your face close to the treatment panel. It is also recommended that you close your eyes during the treatment.

 Press and hold to turn on the Light Therapy Panel. Press again to adjust the Intensity, there are five levels. Press and hold to turn Off.

 Press Treatment Programme to select one of the three programmes:

P1 Skin Rejuvenation - Red Light (625nm)

P2 Deep Skin Rejuvenation - Infrared Light (850nm)

P3 Combined Skin Rejuvenation & Skin Balance - Red & Infrared Light (625nm + 850nm) and Yellow Light (590nm)

 Press to select the required duration (5, 10, 15 or 60 minutes).

 While using the mirrors there is the option of using the mirror light with warm, cool or natural light modes. Press repeatedly to change the light mode. Press and hold to brighten or dim. The Treatment Programme and Mirror Light cannot be used at the same time.

AFTERCARE & REACTIONS

Following a treatment, some individuals may notice temporary skin responses:

- Slight redness can appear post-use but typically fades within a few hours.
- In some cases, users may experience minor peeling, dryness, or itchiness, which usually resolves in a short time.
- Occasional mild headaches have been reported after sessions.

The device uses light wavelengths generally considered safe for cosmetic use. However, if you experience any unusual discomfort or skin reaction, discontinue use immediately and consult a healthcare professional if necessary.

Use of optional treatment boosting serums and gels:

Only use approved or recommended products. A list of approved products is available.
Please contact: help@riobeauty.com.

CLEANING

After each use:

- Power off the device and unplug from any charging source.
- Wipe the Light Therapy Panel with a soft, damp cloth.
- Avoid harsh chemicals, abrasive materials, or submerging the unit in water.

STORAGE

- Store the Light Therapy Panel in a clean, dry place at room temperature.
- Avoid exposure to direct sunlight, excessive humidity, or dust.


DISPOSAL

 The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly

SPECIFICATION

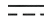Input: 5V  1.0A 5.0W 

Battery: Li-Ion 3.7V 3000mAh 5Wh 11.1Wh

Charging time: 3-3.5 hours

Maximum run time: 2 hours

Model number: MMLT

Program time: 5-60 minutes

Wavelengths:

- Red - 625nm
- Yellow - 590nm
- Infrared - 850nm

This product complies with current EU product legislation.

Read these instructions carefully before use and then retain for future use.

This product is for indoor use only.

Direct current.

This symbol indicates the USB connector.

This unit complies with current UKCA product legislation.

Recycle Paper.

CUSTOMER CARE

If you have any questions, please contact us at: help@riobeauty.com

To reorder additional accessories, please visit www.riobeauty.com or contact your local distributor.

GUARANTEE

To register your two year guarantee online, visit:

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©The Dezac Group Ltd 2026, PO Box 17, Cheltenham Spa, England GL53 7ET
UK RP: The Dezac Group Ltd, PO Box 17,
RP UE: The Dezac Group (Europe) B.V. Azura,
UK Product RP: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, England, GL53 7ET
EU Product RP: The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Netherlands