STYLPRO™

Dual Microcurrent Face & Body Toning Device

Appareil de tonification du visage et du corps par double microcourant
Gesichts- und Körperstraffungs-Gerät mit dualen Mikrostrom
Toningapparaat voor gezicht en lichaam met dubbele microcurrent
Urządzenie do tonizacji twarzy i ciała Dual Microcurrent
Dispositivo tonificador para cuerpo y cara con microcorriente dual
Dispositivo con doppia microcorrente per tonificare viso e corpo
Ansikts- og kroppstoningsapparat med mikrostrømteknologi
Aparelho duplo de microcorrente para tonificação facial e corporal
Kahta mikrovirtaa käyttävä kiinteytyslaite kasvoille ja vartalolle
Apparat med dobbelt mikrostrøm til ansigts- og kropstræning
Ansikts- och kroppstoningsapparat med dubbla mikroströmmar

(EN) USE & CARE GUIDE

FOR YOUR SAFETY AND CONTINUED ENJOYMENT OF THIS PRODUCT, ALWAYS READ INSTRUCTION LEAFLET CAREFULLY BEFORE USING.

The product is intended for private use only and is not intended for commercial use or for use in other applications. This is not a medical product or a medical device.


SPECIFICATIONS:

Name: STYLPRO™ Dual Microcurrent Face & Body Toning Device
Model: MCB01A
Device Size: 60 x 25 x 57 mm
Materials: Abs Plastic, Silicone
Nominal Voltage: 3.7V
Rated capacity of the battery: 500mAh
Energy Capacity of the battery: 1.85Wh
Charge Voltage: 5V
Input power supply: USB socket, 5V 1A
Charging Time: 3 hours
Uses: 30 Uses

INCLUDES:

1x STYLPRO™ Dual Microcurrent Device
1x USB-C Cable
1x Protective Travel Bag

3.7 V  500mAh, 1.85Wh

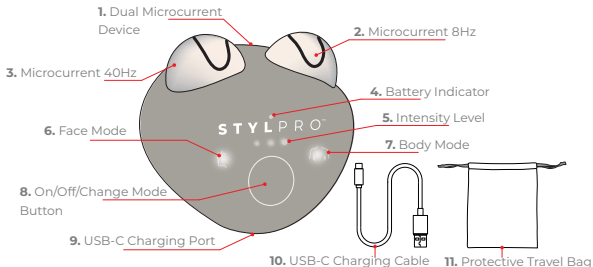MCB01A
Made in China.

©2026 STYLIDEAS Ltd. All rights reserved.
STYLIDEAS owns trademarks, registered trademarks, patents, patents pending, design registrations and protection of the copyright of its products. www.styleideas.com/ip.
Distributed by STYLIDEAS Ltd.
UK: Amshold House, IG10 2RW, UK. +44(0)20 3225 5559 teamtom@styleideas.com.
EU RP: OBELIS S.A. Bd. Général Wahis 53, B-1030 Brussels, Belgium

TEST YOURSELF ON YOUR ARM PRIOR TO USING THE DEVICE ON YOUR FACE AND BODY:

1. PRESS THE ON/OFF BUTTON.
2. PLACE THE DEVICE ON YOUR ARM AND CONDUCT A 3 MINUTE SESSION.
3. AFTER 3 MINUTES, TURN THE DEVICE OFF. VISUALLY CHECK THE AREA. THE AREA MAY INITIALLY TURN RED DUE TO THE HEATING OF YOUR SKIN, BUT IF THE REDNESS LASTS FOR MORE THAN 3 HOURS, YOU MAY HAVE SKIN SENSITIVITY TO MICROCURRENT: SEEK MEDICAL ADVICE BEFORE DECIDING HOW BEST TO USE THIS DEVICE ON YOUR SKIN.

SAFETY INSTRUCTIONS:

1. THIS APPLIANCE IS DESIGNED FOR INDOOR HOME USE ONLY.
2. CONSULT A DOCTOR BEFORE USE, IF YOU HAVE ANY OF THE CONDITIONS BELOW, USE OF THE BODY & FACE SMOOTHER IS NOT RECOMMENDED WITHOUT YOUR DOCTOR'S APPROVAL: METAL IMPLANTS, RECENT SURGERY, OR HAVE A MEDICAL IMPLANT SUCH AS A PACEMAKER, CURRENTLY PREGNANT OR BREASTFEEDING, NAUSEA, HEADACHE, FEVER, SKIN DISEASES, BURNS OR REDDENING, BROKEN BONES, CANCER OR OTHER MALIGNANT CONDITIONS, OSTEOPOROSIS AND MENTAL ILLNESS.
3. FOR ADULT USE ONLY. DO NOT USE THE STYLPRO™ DUAL MICROCURRENT WITH AN INFANT, CHILD, PET, INCAPACITATED OR IMMOBILE PERSON, OR PERSONS WITH REDUCED PHYSICAL OR MENTAL CAPABILITIES WITHOUT SUPERVISION. THIS IS NOT A TOY.
4. IF THE STYLPRO™ DUAL MICROCURRENT MALFUNCTIONS OR IS DAMAGED IN ANY WAY, STOP USING IT IMMEDIATELY. CONTACT STYLPRO™ CUSTOMER SERVICE (SEE SECTION LIMITED TWO-YEAR WARRANTY FOR CONTACT DETAILS).
5. DO NOT USE ON BROKEN, IRRITATED, SWOLLEN, SUNBURNT, OR NUMB/SENSITIVE SKIN, AREAS OF SWELLING, BRUISING, INFLAMMATION, OR SKIN CONDITIONS.
6. DO NOT USE ON THE CHEST / BREAST AREA, EYE AREA, THYROID AREA MID-LINE OF NECK, OR THE GENITALS / GROIN AREA .
7. AVOID THE THYROID AREA: THE MIDDLE AND LOWER PART OF THE CENTER OF YOUR NECK WHERE THE THYROID GLAND IS LOCATED. OVERSTIMULATION OF THE THYROID CAN DISRUPT HORMONE FUNCTIONS.
8. TEST YOURSELF FOR MICROCURRENT SENSITIVITY ON YOUR ARM PRIOR TO USING THE DUAL MICROCURRENT.
9. USE WITH CARE AROUND THE TEETH, EYES, NECK, AND TEMPLES. THESE AREAS HAVE SENSITIVE NERVES. STOP USING IN THESE AREAS IF YOU FEEL PAIN, ODD SENSATIONS, OR FLASHING VISION.

HOW TO USE:

1. ENSURE YOU FACE AND BODY IS MAKEUP FREE AND CLEAN.
2. APPLY YOUR PREFERRED SERUM OR MOISTURISER TO YOUR FACE OR BODY. THIS IS ESSENTIAL TO MAKE SURE THE MICROCURRENTS CONDUCT TO YOUR SKIN. FOR BEST RESULTS, USE A SERUM OR MOISTURISER THAT STAYS WET AND IS NEAR PH NEUTRAL.
3. LONG PRESS THE ON/OFF/MODE BUTTON TO TURN ON THE DEVICE IT WILL START ON FACE MODE.
 - A. TO SWITCH TO BODY MODE, LONG PRESS THE ON/OFF/MODE BUTTON. AN INDICATOR LIGHT SHOULD CHANGE BETWEEN THE FACE OR BODY ICON.
 - B. TO CHANGE THE INTENSITY OF THE MICROCURRENTS, SHORT PRESS THE ON/OFF/MODE BUTTON. AN INDICATOR WILL DISPLAY BETWEEN 1 AND 3 DOTS DENOTING THE INTENSITY LEVELS.
4. LIGHTLY PRESS BOTH METAL SPHERES AGAINST THE SKIN. THE CENTRE OF THE SPHERES PROVIDES A FREQUENCY OF 8 HZ WHILE THE OUTSIDE OF THE SPHERES PROVIDES 40 HZ.
5. GENTLY GLIDE THE DEVICE SLOWLY ACROSS THE SKIN FOR BETWEEN 20 SECONDS AND 3 MINUTES PER AREA, THE MICROCURRENTS SHOULD TINGLE AND TWITCH THE MUSCLES GENTLY BUT IF YOU FIND THEM PAINFUL THEN YOU SHOULD LOWER THE INTENSITY. TOO MUCH INTENSITY MAY GIVE WORSE RESULTS.
 - A. USE ON THE FACE AROUND THE JAWLINE, CHEEKBONES, AND THE FOREHEAD. USE WITH CARE AROUND SENSITIVE AREAS SUCH AS THE TEETH, EYE AREA AND TEMPLES.
 - B. USE ON THE ARMS, LEGS, AND TORSO. AVOID USING NEAR THE HEART.
- 2.

6. LONG PRESS THE BUTTON TO TURN OFF THE DEVICE. WIPE CLEAN WITH A DAMP CLOTH.

FACE MODE IS INDICATED BY AN LED IN THE FACE SYMBOL. THIS MODE PROVIDES TWO SEPARATE MICROCURRENTS. THE INTENSITIES HAVE BEEN OPTIMISED FOR THE FACE AND NECK. CHOOSE FROM 3 MICROCURRENT INTENSITIES.

BODY MODE IS INDICATED BY AN LED IN THE BODY SYMBOL. THIS MODE ALSO PROVIDES TWO SEPARATE MICROCURRENTS. THE INTENSITIES HAVE BEEN OPTIMISED FOR THE BODY. CHOOSE FROM 3 MICROCURRENT INTENSITIES.

DEVICE MOVEMENTS:

WHEN USING THE DEVICE ON YOUR FACE, WE RECOMMEND THESE MOVEMENTS:

1. JAWLINE: FROM CHIN TO OUTER CHEEK. (USE WITH CARE AROUND THE SENSITIVE UPPER CORNERS OF THE MOUTH)
2. CROW'S FEET: FROM CHEEKBONE TO TEMPLES.
3. FOREHEAD: FROM BROWS TO HAIRLINE. (USE WITH CARE AROUND THE SENSITIVE TEMPLES)
4. NECK: FROM SLIGHTLY OFF-CENTRE, AVOIDING THE MIDLINE OF THE NECK. GLIDE THE DEVICE TOWARDS THE BACK OF THE NECK.

WHEN USING THE DEVICE ON YOUR BODY, WE RECOMMEND THESE MOVEMENTS:

1. FOREARMS: FROM WRIST TO ELBOW.
2. UPPER ARMS: FROM ELBOW TO SHOULDER.
3. SHOULDERS: FROM UPPER TRAPEZIUS OUT TOWARD SHOULDER TIP.
4. SHOULDERS: FROM COLLARBONE SWEEPING UP AND OVER THE SHOULDER.
5. TORSO: FROM LOWER ABDOMEN UP TOWARD THE RIBS (LEFT AND RIGHT SIDES).
6. TORSO: FROM HIPS UP TO WAIST, THEN WAIST UP TO UNDER-BUST/RIB LINE.
7. SIDES: FROM OUTER HIP, UP TOWARD THE UNDERARM.
8. LEGS: FROM ANKLE UP TO KNEE, THEN KNEE UP TO THIGH AND TOWARD THE HIP/GLUTE.
9. GLUTES: FROM THE BASE OF THE GLUTE TO THE TOP OF THE GLUTE.

HOW TO CLEAN:

1. USE A DISINFECTANT WIPE TO WIPE DOWN THE DEVICE AND REMOVE ANY SERUM. WE RECOMMEND TO CLEAN AFTER EVERY USE.

CHARGING:

1. INSERT THE USB-C CABLE (PROVIDED) INTO THE PORT AT THE BOTTOM OF THE DEVICE.
2. ATTACH TO A USB FRIENDLY POWER SOURCE.
3. THE WHITE LEDS ON THE FRONT OF THE DEVICE WILL FLASH WHEN CHARGING.
4. ONCE FULLY CHARGED A GREEN LED WILL SHINE AND THE WHITE LEDS WILL STOP FLASHING.

SAFETY CONSIDERATIONS:

Use of a Microcurrent skincare device in any way other than those described in this manual voids the warranty and could result in serious injury.

The STYLPRO™ DUAL MICROCURRENT is an electrical device. To avoid electric shock, adhere to the following safety considerations:

- Do not get your device wet or submerge your device in water.
- Do not wash electrical parts with water or liquids.
- Do not directly touch the device or switches with wet hands.

- Do not use the device if it has been dropped, left outside, or has encountered water or liquids.
- Do not modify, or attempt to modify or repair, this equipment. There are no user serviceable parts.
- Do not use or store your device in hot temperatures that are or can exceed 54°C (130°F).

WARNING:

TO REDUCE THE RISK OF BURNS, ELECTROCUTION, FIRE, OR INJURY TO PERSONS:

1. READ ALL INSTRUCTIONS BEFORE USING THE APPLIANCE.
2. IF YOU FEEL ABNORMALITIES IN YOUR SKIN DURING USE OR AFTER USE, PLEASE CONSULT YOUR MEDICAL SPECIALIST.
3. TO PROTECT AGAINST ELECTRICAL HAZARDS, DO NOT USE WHILE BATHING OR IN THE SHOWER.
4. DO NOT OPERATE THIS APPLIANCE IF THE APPLIANCE MALFUNCTIONS OR HAS BEEN DROPPED OR DAMAGED IN ANY MANNER. RETURN APPLIANCE TO NEAREST SERVICE FACILITY FOR EXAMINATION, ELECTRICAL OR MECHANICAL ADJUSTMENT, OR REPAIR. NO REPAIRS SHOULD BE ATTEMPTED BY THE CONSUMER.
5. USE APPLIANCE ONLY FOR INTENDED USE AS DESCRIBED IN THIS INSTRUCTION BOOK. THE USE OF ATTACHMENTS NOT RECOMMENDED BY STYLPRO™ MAY CAUSE HAZARDS.
6. DO NOT PLACE ANY PART OF THE PRODUCT NEAR HEAT SOURCES SUCH AS RADIATORS, COOKERS, ETC DO NOT LAY ON A SOFT SURFACE SUCH AS A BED AS DEVICE MAY OVERHEAT.
7. KEEP THE DEVICE AWAY FROM CORROSIVE LIQUIDS.
8. DO NOT SHARE THIS PRODUCT WITH ANYONE ELSE.
9. PLEASE DO NOT USE DURING CHARGING.
10. DO NOT PLACE OR STORE WHERE IT CAN FALL OR BE PULLED INTO A TUB OR SINK. IF IT FALLS INTO WATER, UNPLUG IMMEDIATELY—DO NOT REACH INTO WATER
11. DO NOT USE DIRECTLY OVER OR NEAR METAL IMPLANTS, DENTAL FILLINGS, BRACES, OR OTHER METALLIC DENTAL WORK, AS THIS MAY CAUSE DISCOMFORT, PAIN, OR A METALLIC TASTE.

BATTERIES: THE PRODUCT CONTAINS LITHIUM-ION BATTERY. RISK OF FIRE AND BURNS, DO NOT OPEN, CRUSH, HEAT ABOVE 60°C/140 °F OR INCINERATE. DO NOT SHORT CIRCUIT. IF THE LITHIUM-ION BATTERY BULGES SEVERELY, DISCONTINUE USE. KEEP SMALL CELLS WHICH ARE CONSIDERED SWALLOWABLE OUT OF REACH OF CHILDREN. SWALLOWING MAY LEAD TO BURNS, PERFORATION OF SOFT TISSUE, AND DEATH. SEVERE BURNS CAN OCCUR WITHIN 2 HOURS OF INGESTION. IN CASE OF INGESTION OF A CELL, SEEK MEDICAL ATTENTION PROMPTLY. NEVER RECHARGE NON-RECHARGEABLE BATTERIES. THIS PRODUCT HAS A BUILT-IN RECHARGEABLE BATTERY WHICH CANNOT BE REPLACED BY THE USER. THE REMOVAL OR REPLACEMENT OF THE RECHARGEABLE BATTERY MAY ONLY BE CARRIED OUT BY THE MANUFACTURER OR HIS CUSTOMER SERVICE OR BY A SIMILARLY QUALIFIED PERSON IN ORDER TO AVOID HAZARDS. WHEN DISPOSING OF THE PRODUCT, IT SHOULD BE NOTED THAT THIS PRODUCT CONTAINS A RECHARGEABLE BATTERY.

STORAGE:

ENSURE DEVICE IS 'OFF' BEFORE STORAGE. ENSURE THE DEVICE IS UNPLUGGED FROM ANY POWER SUPPLY. PLACE IN COOL, DRY AREA OUT OF REACH OF CHILDREN, WHERE IT WILL NOT BE DAMAGED.

MAINTENANCE:

YOUR STYLPRO™ DUAL MICROCURRENT IS VIRTUALLY MAINTENANCE FREE. TO CLEAN TURN THE APPLIANCE OFF, LET THE DEVICE COOL TO ROOM TEMPERATURE, AND WIPE CLEAN WITH A DAMP CLOTH. IF ANY ABNORMAL CONDITION OCCURS, TURN OFF & CONTACT STYLPRO™ CUSTOMER SERVICES ON +44 (0) 20 3225 5559.

LIMITED TWO-YEAR WARRANTY:

OUR GOODS COME WITH GUARANTEES THAT CANNOT BE EXCLUDED UNDER THE EUROPEAN CONSUMER LAW. EU DIRECTIVE 1999/44/EC, THIS DIRECTIVE STATES 'A TWO-YEAR GUARANTEE APPLIES FOR THE SALE OF ALL

4.

CONSUMER GOODS EVERYWHERE IN THE EU. YOU ARE ENTITLED TO A REPLACEMENT OR REFUND FOR A MAJOR FAILURE YOU ARE ALSO ENTITLED TO HAVE THE GOODS REPAIRED OR REPLACED IF THE GOODS FAIL TO BE OF ACCEPTABLE QUALITY AND THE FAILURE DOES NOT AMOUNT TO A MAJOR FAILURE. STYLPRO™ WILL REPAIR OR REPLACE (AT OUR OPTION) YOUR APPLIANCE FREE OF CHARGE FOR TWENTY-FOUR MONTHS FROM THE DATE OF PURCHASE, IF THE APPLIANCE IS DEFECTIVE IN WORKMANSHIP OR MATERIALS. THIS WARRANTY IS AN ADDITIONAL BENEFIT AND DOES NOT AFFECT YOUR LEGAL RIGHTS. TO OBTAIN SERVICE UNDER THIS WARRANTY, PLEASE CONTACT STYLPRO™ CUSTOMER SERVICE AT TEAMTOM@STYLIDEAS.COM +44 (0) 20 3225 5559. TO MAKE A WARRANTY CLAIM YOU MUST RETAIN YOUR PROOF OF PURCHASE.

IMPORTANT: THIS WARRANTY DOES NOT COVER DAMAGE RESULTING FROM ACCIDENTS, MISUSE, ALTERATIONS, ABUSE, LACK OF REASONABLE CARE OR USE OTHER THAN DESCRIBED IN THESE INSTRUCTIONS, OR PERFORMANCE OF SERVICE OR REPAIRS BY AN UNAUTHORISED SERVICE CENTRE. THIS WARRANTY IS VALID ONLY ON MERCHANDISE SOLD BY AUTHORISED RETAILERS. NOTE, THE ITEM MUST BE THOROUGHLY CLEANED PRIOR TO RETURN.

CORRECT DISPOSAL OF THIS PRODUCT:

Disposal of Product: The WEEE symbol on the STYLIDEAS, the Original Accessories or packaging indicates that this device must not be treated as unsorted municipal waste but must be collected separately! Dispose of the STYLIDEAS Product and its Original Accessories via a collection point for the recycling of waste electrical and electronic equipment if you live within the EU and in other European countries that operate separate collection systems for waste electrical and electronic equipment. By disposing of the STYLIDEAS Product and its Original Accessories in the proper manner, you help to avoid possible hazards for the environment and public health that could otherwise be caused by improper treatment of waste equipment. The recycling of materials contributes to the conservation of natural resources. Therefore, do not dispose of your old electrical and electronic equipment with the unsorted municipal waste.

Disposal of Packaging Waste:

The packaging is made of environmentally friendly materials, which may be disposed through your local recycling facilities. By disposing of the packaging in the proper manner, you help to avoid possible hazards for the environment and public health. The symbol on the carton packaging indicates that the packaging is made of PAP.

Disposal of Batteries:

The STYLIDEAS Product contains a battery. Batteries may not be disposed of with the usual domestic waste. They may contain toxic heavy metals and are subject to hazard our waste regulations. For this reason, dispose the Rechargeable STYLIDEAS Product (the part that contains the battery) at a local collection point or the recycling of waste electrical and electronic equipment. Do not attempt to open the STYLIDEAS Product and to remove the battery yourself!

CE MARKING: UK/EU – DECLARATION OF CONFORMITY: STYLIDEAS LTD DECLARES THAT STYLPRO™ DUAL MICROCURRENT COMPLIES WITH THE ESSENTIAL REQUIREMENTS AND OTHER RELEVANT PROVISIONS. A COPY OF THE DECLARATION OF CONFORMITY IS AVAILABLE ON REQUEST. STYLIDEAS LTD, AMSHOLD HOUSE, GOLDINGS HILL, LOUGHTON, ESSEX, IG10 2RW, UK. OBELIS S.A, BD. GÉNÉRAL WAHIS 53, B-1030 BRUSSELS, BELGIUM.

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- **CE-MÄRKNING: UK/EU – FÖRSÄKRAN OM ÖVERRENSSTÄMMELSE:** STYLIDEAS LTD FÖRKLARAR ATT STYLPRO™
- Kopia av FÖRSÄKRAN OM ÖVERRENSSTÄMMELSE ÄR TILLGÄNGLIG VID BEGÄRAN. STYLIDEAS LTD, AMSHOLD
- HOUSE, GOLDINGS HILL, LOUGHTON, ESSEX, IG10 2RW, STORBRITANNIEN. OBELIS S.A, BD. GÉNÉRAL WAHIS 53,