SCAN ME

©2026 STYLIDEAS Ltd. All rights reserved. STYLIDEAS owns trademarks, registered trademarks, patents, patents pending, design registrations and protection of the copyright of its products. www.stylideas.com/ip.

Distributed by STYLIDEAS Ltd,

UK: Amshold House, IG10 2RW, UK. +44(0)20 3225 5559
teamtom@stylideas.com.

EU: OBELIS S.A, Bd. Général Wahis 53, B-1030 Brussels, Belgium

For indoor use only

Li-ion

RoHS COMPLIANT

PAP

Made in China.

RS01A 56

STYLPRO™

SPEC-TACULAR

EMS & Red LED Under Eye Glasses

Lunettes pour zone sous les yeux à EMS et LED rouge ■ Gafas de luz LED roja y EMS para la zona inferior del ojo ■ EMS- & Rotlicht-LED-Brille unter den Augen ■ Bril met EMS en rood licht ■ Okulary pod oczy z EMS i czerwonymi diodami LED ■ Occhiali con EMS & LED a luce rossa per la cura della pelle sotto gli occhi ■ EMS- og LED-briller for rødlysterapi ■ Óculos para a zona debaixo dos olhos com EMS e LED vermelho ■ EMS- og LED-rødlysbriller ■ EMS- ja LED-punavalolasis ■ Glasögon med EMS och röd LED för området under ögat

(EN) USE & CARE GUIDE

FOR YOUR SAFETY AND CONTINUED ENJOYMENT OF THIS PRODUCT, ALWAYS READ INSTRUCTION LEAFLET CAREFULLY BEFORE USING.

The product is intended for private use only and is not intended for commercial use or for use in other applications. This is not a medical product or a medical device.

SPECIFICATIONS

Name: STYLPRO™ Spec-Tacular EMS & Red LED Under Eye Glasses

Model: RS01A

Device Size: Folded Size 140mm x 55mm x 53mm.

Materials: ABS, Rubber & Silicone

Nominal Voltage: 3.7V

Rated capacity of the battery: 90mAh

Energy Capacity of the battery: 0.33Wh

Charge Voltage: 4.2V


Input power supply: USB socket, 5V, 1A

Charging Time: 1.5 hours

3.7 V , 90mAh, 0.33Wh

PARTS IDENTIFICATION:

1. Red Light Spectacles
2. On/Off Button
3. Intensity Button
4. USB Charging Point
5. Red Light Pads

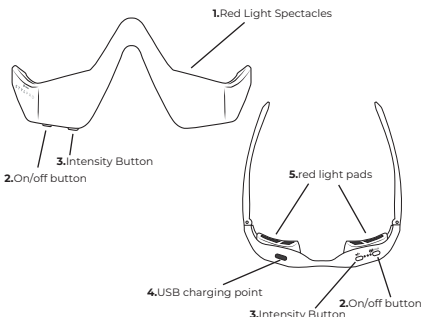SAFETY GUIDELINES:

1. Please do not use if you have eye inflammation or conjunctivitis.
2. Only for use on the eye area.
3. We do not advise using this device if you are epileptic, photophobic, or light sensitive, if you are pregnant or have had skin cancer in the last 12 months, have undergone facial plastic surgery, have a metal implant, are pregnant, or if you are currently on your menstrual cycle.
4. Do not use this product when wearing contact lenses.
5. **DO NOT USE THE PRODUCT WITHOUT EYE CREAM OR SERUM.**

6. For adult use only. Not intended for use by children. This is not a toy.
7. Please wipe off any excess skin care products which may accumulate, to clean, wipe gently with a damp cloth or cleaning cloth provided.
8. Avoid storing the STYLPRO™ Spec-tacular Glasses in hot, humid, or damp areas.
9. If the STYLPRO™ Spec-tacular Glasses malfunctions or is damaged in any way, stop using it immediately. Contact STYLPRO™ CUSTOMER SERVICE (See section LIMITED TWO-YEAR WARRANTY for contact details).
10. Use LED lights with care. Extended use and prolonged exposure of LED lights into the eyes may cause damage. LEDs are not replaceable.
11. Do not use on sensitive, irritated or broken skin. Do not use on swollen or sunburnt skin or if you have any other skin conditions.
12. Do Not look directly into the LED lights.

CHARGING:

1. Insert the USB cable (provided) into the port at the bottom of the glasses. Attach to a USB friendly power source. While charging, the indicator lights will flash red. When the device is fully charged the indicator lights will turn green.

HOW TO USE:

Removal of any make-up is advisable before use.

1. Remove your makeup with an oil-free makeup remover or cleanser.
2. Apply your favourite serum to your under eye area. EMS is less effective on dry skin. Do not use the device without serum or if the serum has dried out on the skin.
3. Long press the On/Off button to turn on the device. You will see 1 red LED light turn on at the bottom of the device.
4. Short press the On/Off button to switch between EMS frequencies. There are 2 frequency settings, continuous and pulse.
5. To adjust the EMS intensity, short press the intensity button. There are 3 intensity settings. **ALWAYS START ON INTENSITY SETTING 1.**
6. Place the device on to your under eye area and relax. (If you do not feel any sensation at all, apply more serum to the skin and or press the device to the skin.)
7. Long press the On/Off button to turn off the device. Wipe the device clean with a slightly damp cloth.

THE DEVICE WILL AUTOMATICALLY TURN OFF AFTER A 10 MINUTE PERIOD.

IT IS NORMAL TO FEEL LIGHT ELECTRIC SHOCK SENSATION WHEN USING THE DEVICE. THE MORE MOIST THE SKIN, THE STRONGER THE FEELING. YOU SHOULD NOT FEEL DISCOMFORT OR PAIN. IF YOU DO THEN LOWER THE SETTING OR ADJUST THE GLASSES TO A MORE COMFORTABLE POSITION.

SAFETY INFORMATION :

Red Light products are not intended to diagnose, treat, cure, or prevent any disease. Our products are low-risk, general wellness / fitness products that do not require medical clearance. They are designed as "General Wellness" home use products only.

SAFETY CONSIDERATIONS:

Use of a Red Light therapy device in any way other than those described in this manual voids the warranty and could result in serious injury.

STYLPRO™ Spec-tacular Glasses are electrical devices. To avoid electric shock, adhere to the following safety considerations:

- Do not get your device wet or submerge your device in water.
- Do not wash electrical parts with water or liquids.
- Do not directly touch the device or switches with wet hands.

- Do not pull, carry, or lift the device by its charging cord.
- If the charging cord is damaged, do not use the device.
- Do not use the device if it has been dropped, left outside, or has encountered water or liquids.
- Do not use on or near heated surfaces.
- Do not modify, or attempt to modify or repair, this equipment. There are no user serviceable parts.
- Do not use or store your device in hot temperatures that are or can exceed 54°C (130°F).

IMPORTANT SAFETY INSTRUCTIONS: WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE TAKEN, INCLUDING THE FOLLOWING: READ ALL INSTRUCTIONS BEFORE USING APPLIANCE. SAVE THESE INSTRUCTIONS. KEEP AWAY FROM WATER. THIS PRODUCT IS NOT A TOY. DO NOT STARE DIRECTLY INTO THE LIGHTS. USE INDOORS ONLY. FOR DOMESTIC HOUSEHOLD USE ONLY. DO NOT ATTEMPT TO OPEN THE BACK OF THE DEVICE. IF REQUIRED, CONSULT A PROFESSIONAL. KEEP AWAY FROM SOURCES OF DIRECT HEAT AND SUNLIGHT. KEEP AWAY FROM CHILDREN.

DANGER: TO REDUCE THE RISK OF ELECTRIC SHOCK: ALWAYS TURN OFF BEFORE CLEANING OR STORING.

WARNING:

TO REDUCE THE RISK OF BURNS, ELECTROCUTION, FIRE, OR INJURY TO PERSONS:

1. READ ALL INSTRUCTIONS BEFORE USING THE APPLIANCE
2. THIS PRODUCT SHOULD NOT BE USED BY; CHILDREN, PREGNANT WOMEN, PATIENTS WITH EPILEPSY, TUMORS AND CANCER PATIENTS, PATIENTS WITH CARDIAC PACEMAKERS OR SIMILAR.
3. REMOVE , GLASSES, FALSE EYELASHES AND CONTACT LENSES BEFORE USE.
4. IF YOU FEEL ABNORMALITIES IN YOUR SKIN DURING USE OR AFTER USE, PLEASE CONSULT YOUR MEDICAL SPECIALIST.
5. TO PROTECT AGAINST ELECTRICAL HAZARDS, DO NOT USE WHILE BATHING OR IN THE SHOWER.
6. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED BY OR NEAR CHILDREN OR PERSONS MADE WEAK OR DISABLED BY ILLNESS OR INJURY.
7. DO NOT USE ON SENSITIVE OR BROKEN SKIN. DO NOT USE ON SWOLLEN OR SUNBURNT SKIN OR IF YOU HAVE ANY OTHER SKIN CONDITIONS.
8. DO NOT OPERATE THIS APPLIANCE IF THE APPLIANCE MALFUNCTIONS OR HAS BEEN DROPPED OR DAMAGED IN ANY MANNER. RETURN APPLIANCE TO NEAREST SERVICE FACILITY FOR EXAMINATION, ELECTRICAL OR MECHANICAL ADJUSTMENT, OR REPAIR. NO REPAIRS SHOULD BE ATTEMPTED BY THE CONSUMER.
9. DO NOT DIRECTLY TOUCH THE EYEBALL.
10. USE APPLIANCE ONLY FOR INTENDED USE AS DESCRIBED IN THIS INSTRUCTION BOOK. THE USE OF ATTACHMENTS NOT RECOMMENDED BY STYLPRO™ MAY CAUSE HAZARDS.
11. DO NOT USE OUTDOORS OR OPERATE WHERE AEROSOL (SPRAY) PRODUCTS ARE BEING USED, OR WHERE OXYGEN IS BEING ADMINISTERED.
12. DO NOT PLACE ANY PART OF THE PRODUCT NEAR HEAT SOURCES SUCH AS RADIATORS, COOKERS, ETC DO NOT LAY ON A SOFT SURFACE SUCH AS A BED AS DEVICE MAY OVERHEAT.
13. KEEP THE DEVICE AWAY FROM CORROSIVE LIQUIDS
14. ENSURE DEVICE IS 'OFF' BEFORE STORAGE.
15. DO NOT USE APPLIANCE IN THE BATHROOM OR NEAR WATER.
16. DO NOT PLACE OR STORE APPLIANCE WHERE IT CAN FALL OR BE PULLED INTO A TUB OR SINK.
17. DO NOT PLACE IN OR DROP INTO WATER OR OTHER LIQUIDS.

18. IF THE APPLIANCE FALLS INTO WATER, DO NOT REACH INTO THE WATER.
19. CHOKING HAZARD DUE TO SMALL PARTS. NOT SUITABLE FOR CHILDREN UNDER 3 YEARS.
20. DO NOT SHARE THIS PRODUCT WITH ANYONE ELSE.
21. PLEASE DO NOT USE DURING CHARGING.

BATTERIES: THE PRODUCT CONTAINS LITHIUM-ION BATTERY. RISK OF FIRE AND BURNS, DO NOT OPEN, CRUSH, HEAT ABOVE 60°C/140 °F OR INCINERATE. DO NOT SHORT CIRCUIT. IF THE LITHIUM-ION BATTERY BULGES SEVERELY, DISCONTINUE USE. KEEP SMALL CELLS WHICH ARE CONSIDERED SWALLOWABLE OUT OF REACH OF CHILDREN. SWALLOWING MAY LEAD TO BURNS, PERFORATION OF SOFT TISSUE, AND DEATH. SEVERE BURNS CAN OCCUR WITHIN 2 HOURS OF INGESTION. IN CASE OF INGESTION OF A CELL, SEEK MEDICAL ATTENTION PROMPTLY. NEVER RECHARGE NON-RECHARGEABLE BATTERIES. THIS PRODUCT HAS A BUILT-IN RECHARGEABLE BATTERY WHICH CANNOT BE REPLACED BY THE USER. THE REMOVAL OR REPLACEMENT OF THE RECHARGEABLE BATTERY MAY ONLY BE CARRIED OUT BY THE MANUFACTURER OR HIS CUSTOMER SERVICE OR BY A SIMILARLY QUALIFIED PERSON IN ORDER TO AVOID HAZARDS. WHEN DISPOSING OF THE PRODUCT, IT SHOULD BE NOTED THAT THIS PRODUCT CONTAINS A RECHARGEABLE BATTERY.

STORAGE

ENSURE DEVICE IS 'OFF' BEFORE STORAGE. ENSURE THE GLASSES ARE UNPLUGGED FROM ANY POWER SUPPLY. PLACE IN COOL, DRY AREA OUT OF REACH OF CHILDREN, WHERE GLASSES WILL NOT BE DAMAGED.

MAINTENANCE

PLEASE WIPE OFF ANY EXCESS SKIN CARE PRODUCTS WHICH MAY ACCUMULATE.

YOUR STYLPRO™ SPEC-TACULAR GLASSES ARE VIRTUALLY MAINTENANCE FREE. TO CLEAN TURN THE APPLIANCE OFF AND WIPE EXTERIOR WITH A DAMP CLOTH OR CLEANING CLOTH PROVIDED. IF ANY ABNORMAL CONDITION OCCURS, TURN OFF & CONTACT STYLPRO™ CUSTOMER SERVICES ON +44 (0) 20 3225 5559.

LIMITED TWO-YEAR WARRANTY

OUR GOODS COME WITH GUARANTEES THAT CANNOT BE EXCLUDED UNDER THE EUROPEAN CONSUMER LAW. EU DIRECTIVE 1999/44/EC, THIS DIRECTIVE STATES 'A TWO-YEAR GUARANTEE APPLIES FOR THE SALE OF ALL CONSUMER GOODS EVERYWHERE IN THE EU. YOU ARE ENTITLED TO A REPLACEMENT OR REFUND FOR A MAJOR FAILURE YOU ARE ALSO ENTITLED TO HAVE THE GOODS REPAIRED OR REPLACED IF THE GOODS FAIL TO BE OF ACCEPTABLE QUALITY AND THE FAILURE DOES NOT AMOUNT TO A MAJOR FAILURE. STYLPRO™ WILL REPAIR OR REPLACE (AT OUR OPTION) YOUR APPLIANCE FREE OF CHARGE FOR TWENTY-FOUR MONTHS FROM THE DATE OF PURCHASE, IF THE APPLIANCE IS DEFECTIVE IN WORKMANSHIP OR MATERIALS. THIS WARRANTY IS AN ADDITIONAL BENEFIT AND DOES NOT AFFECT YOUR LEGAL RIGHTS. TO OBTAIN SERVICE UNDER THIS WARRANTY, PLEASE CONTACT STYLPRO™ CUSTOMER SERVICE AT TEAMTOM@STYLIDEAS.COM +44 (0) 20 3225 5559. TO MAKE A WARRANTY CLAIM YOU MUST RETAIN YOUR PROOF OF PURCHASE.

IMPORTANT: THIS WARRANTY DOES NOT COVER DAMAGE RESULTING FROM ACCIDENTS, MISUSE, ALTERATIONS, ABUSE, LACK OF REASONABLE CARE OR USE OTHER THAN DESCRIBED IN THESE INSTRUCTIONS, OR PERFORMANCE OF SERVICE OR REPAIRS BY AN UNAUTHORISED SERVICE CENTRE. THIS WARRANTY IS VALID ONLY ON MERCHANDISE SOLD BY AUTHORISED RETAILERS. NOTE, THE ITEM MUST BE THOROUGHLY CLEANED PRIOR TO RETURN.

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PRODUCT HAS REACHED THE END OF ITS USAGE, PLEASE DO NOT DISPOSE WITH REGULAR HOUSEHOLD WASTE. RECYCLING IS THE SAFEST AND MOST RESPONSIBLE WAY TO DISPOSE OF ELECTRICAL ITEMS THAT COME TO THE END OF THEIR LIFECYCLE. PLEASE DISPOSE OF IT ACCORDING TO YOUR LOCAL REGULATIONS. THE WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE) DIRECTIVE IN EUROPEAN COUNTRIES STIPULATES CORRECT RECYCLING PRACTICES. DO NOT DISPOSE IN REGULAR HOUSEHOLD WASTE.

UKCA & CE MARKING: UK/EU – DECLARATION OF CONFORMITY: STYLIDEAS LTD DECLARES THAT STYLPRO™ SPEC-TACULAR GLASSES COMPLY WITH THE ESSENTIAL REQUIREMENTS AND OTHER RELEVANT PROVISIONS. A COPY OF THE DECLARATION OF CONFORMITY IS AVAILABLE ON REQUEST. STYLIDEAS LTD, AMSHOLD HOUSE, GOLDINGS HILL, LOUGHTON, ESSEX, IG10 2RW, UK. OBELIS S.A, BD. GÉNÉRAL WAHIS 53, B-1030 BRUSSELS, BELGIUM.

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